

User Manual

Version 1.1 - UMFC3301M



Evercode™ Mouse

TCR/BCR Cell Fixation v3

Mid-Throughput Workflow

For use with

ECFC3301, ECFC3303

ECFC3305, ECFC3501

ECFC3503, ECFC3505

Legal Notices

This document and its contents are proprietary to Parse Biosciences, Inc. ("Parse Biosciences") and are intended solely for use by its customers in connection with the use of the product(s) described herein and for no other purpose. The products may be used solely FOR RESEARCH PURPOSES, AND MAY NOT BE USED IN ANY DIAGNOSTIC OR THERAPEUTIC USE IN HUMANS OR ANIMALS. This document and its contents shall not be used or distributed for any other purpose and/ or otherwise communicated, disclosed or reproduced in any way whatsoever without the prior written consent of Parse Biosciences.

No rights are granted under this document with respect to any of Parse Biosciences' intellectual property rights. The license to use of any products described herein is subject to a separate written agreement between Parse Biosciences and the applicable user.

The instructions in this document must be strictly and explicitly followed by qualified and properly trained personnel in order to ensure the proper and safe use of the product(s) described herein. Parse Biosciences shall have no liability for any direct, indirect, consequential or incidental damages arising out of any failure to use the product(s) in strict compliance with the terms herein.

This document may contain references to third-party sources of information, hardware or software, products, or services and/ or third-party web sites (collectively "Third Party Information"). Parse Biosciences does not control, and is not responsible for, any Third Party Information. The inclusion of Third Party Information in this document does not imply endorsement by Parse Biosciences of the Third Party Information or the third party in any way.

The product(s) described in this document are provided for one-time use by the purchaser and may not be re-used, refurbished or resold. In addition, such product(s) may not be altered, changed or modified by anyone other than Parse Biosciences and its authorized agents, and Parse Biosciences will not be liable for any such alterations, changes or modifications. The product(s) described herein may be covered by one or more of the following patents:

U.S. Pat. No. 10,900,065

U.S. Pat. No. 11,168,355

U.S. Pat. No. 11,427,856

U.S. Pat. No. 11,634,751

U.S. Pat. No. 11,639,519

U.S. Pat. No. 11,680,283

Patents pending in the U.S. and other countries

PARSE and EVERCODE are trademarks of Parse Biosciences, Inc. The information made available herein is the sole property of Parse Biosciences and its licensors.

Copyright (c) 2024 Parse Biosciences, Inc. All Rights Reserved.

Table of Contents

Overview	4
Workflow	4
Protocol Timing.....	6
Important Guidelines	7
Parts List.....	12
User Supplied Equipment and Consumables	14
Section 1: Set Up	17
1.1. Block Tubes with BSA.....	17
1.2. Prepare Master Mixes	19
Section 2: Fixation	22
2.1. Cell Fixation	22
Appendix: Revision History	27

Overview

Workflow

From a single cell suspension, the Evercode Mouse TCR/BCR Cell Fixation kit generates fixed and permeabilized cells ready for use in Evercode TCR and Evercode BCR kits.

This mid-throughput plate-based workflow enables more samples to be fixed in parallel, streamlining fixation when processing >12 samples at a time.

The mid-throughput plate-based workflow is compatible with both the 12 and 96 Reactions Cell Fixation kits.

If processing between 12 and 48 samples we suggest using ECFC3301 kit(s).

If processing 48 samples at a time, we suggest using half of the high throughput ECFC3501 kit(s). Please note that each high throughput kit should not be used more than twice.

Fixation maintains cell structure, prevents RNA degradation, and locks the RNA inside the cells, which are crucial for downstream processing with Evercode split-pool combinatorial barcoding technology (Figure 1).

Because fixed samples are also stable for up to 6 months at -80°C, Evercode Mouse TCR/BCR Cell Fixation provides flexibility by separating sample collection from library preparation. It also enables samples to be stored and batched after fixation so they can be processed through library preparation together, reducing the potential of batch effects.

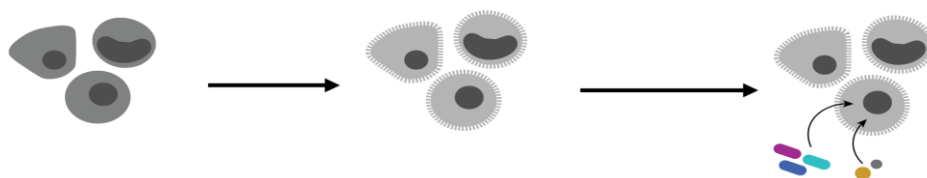


Figure 1: Evercode Mouse TCR/BCR Cell Fixation. Cells in suspension are fixed and permeabilized before undergoing the split-pool combinatorial barcoding steps.

The figure below provides an overview of the fixation workflow. Between 100,000 and 1 million cells can be fixed in a single reaction. Note that more than 100,000 cells may need to be fixed to fully utilize the capacity of the downstream Evercode TCR and BCR kits. See Important Guidelines for additional details.

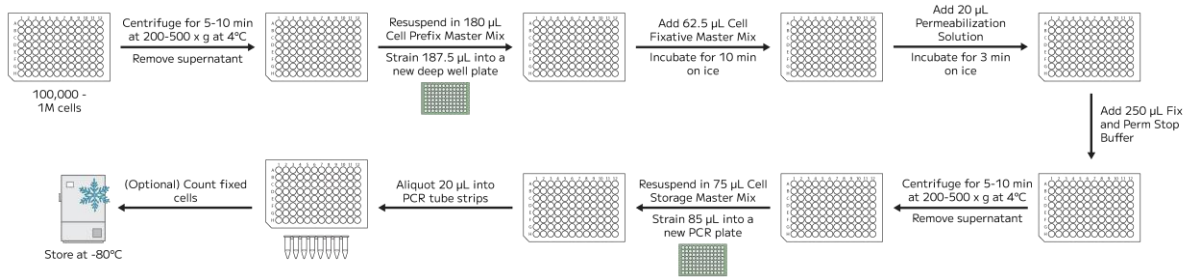


Figure 2: The Mid-Throughput Plate-Based Fixation Evercode Mouse TCR/BCR Cell Fixation Workflow

Protocol Timing

The table below provides details of the total and hands-on time required for the cell fixation workflow.

SECTION	TOTAL TIME	HANDS-ON TIME	STOPPING POINTS
Section 1: Set Up			
1.1 Block plates with BSA	65 min	5 min	4°C ≤ 1 month
1.2 Prepare Master Mixes	15 min	15 min	
Section 2: Fixation			
2.1 Cell Fixation	60 min	60 min	-80°C ≤ 6 months

Important Guidelines

These guidelines provide additional information to obtain optimal performance beyond the detailed instructions in the protocol. For additional questions not discussed below, please contact us at support@parsebiosciences.com. We also have a library of additional resources and videos on our support site at <https://support.parsebiosciences.com/>.

Sample Input

- This protocol begins with a previously prepared single cell suspension. We recommend suspensions with >70% viability (ideally above 90%) and <5% aggregation/debris.
- If cells were previously frozen, ensure the suspension is completely thawed and in suspension before beginning fixation.
- We recommend minimizing the length of time samples are stored on ice prior to fixation, as it can negatively impact results.
- Between 100,000 and 1 million cells can be fixed in a single reaction. However, we recommend using the highest number available up to 1 million total. Exceeding 1 million cells in a single fixation will result in substantially elevated doublet rates.
- The minimum input into fixation should also be determined based on how the samples will be processed downstream. The table below provides guidance on the post-fixation concentrations needed for downstream kits. However, more or less sample input may be required depending on the exact experimental design. To accurately determine required post-fixation cell concentrations and volumes, reference the relevant [Sample Loading Table](#).
- Note that retention during fixation varies typically between 40-60%, and some cells will be lost when freezing and thawing fixed samples, typically between 5-15%. The final concentration of cells post-fixation is also influenced by the resuspension volume used in Step 18 of Section 2. These factors should all be taken into account when determining how much sample input is needed for fixation.

CELL CONCENTRATIONS		
Kit	Target Post-Fixation Concentration	Minimum Post-Thaw Concentration to Fully Load Kit
Evercode TCR or BCR Mini	≥ 500 cells/μL	298 cells/μL
Evercode TCR or BCR	≥ 1,000 cells/μL	520 cells/μL

CELL CONCENTRATIONS		
Evercode TCR or BCR Mega	≥ 3,000 cells/μL	2,126 cells/μL

Avoiding RNase Contamination

- Standard precautions should be taken to avoid introducing RNases into samples or reagents throughout the workflow. Always wear proper laboratory gloves and use aseptic technique.
- Although RNases are not inactivated by ethanol or isopropanol, they are inactivated by products such as RNaseZap RNase Decontamination Solution (Thermo Fisher Scientific). These can be sprayed on benchtops and pipettes.
- Nuclease-free, filtered pipette tips should be used to reduce RNase contamination from pipettes.

Cell Detachment

- If using adherent cell line samples, we recommend TrypLE Express Enzyme (1X), phenol red (Thermo Fisher Scientific). Due to high RNase activity, we do not recommend dissociation with standard trypsin, which may reduce gene and transcript detection.

Cell Counting and Quality Assessment

- We recommend a hemocytometer for counting, but alternative counting devices can also be used. If possible, validate counts from alternative devices to a hemocytometer when first using an Evercode Fixation workflow.
- When first using an Evercode Fixation workflow, we suggest saving images at each counting step.
- To assess sample quality, we recommend using viability stains like trypan blue or acridine orange and propidium iodide (AO/PI).
- After fixation, the cells are permeabilized and should appear dead with viability stains. If using Acridine Orange/Propidium Iodide (AO/PI) stains, we suggest using the red (PI) channel to count to avoid the impact of any autofluorescence in the green (AO) channel.
- Examples of trypan blue stained fixed cells are shown below. High quality fixed samples have single distinct cells with <5% cell aggregation and no debris. Higher levels of aggregation will lead to elevated doublets after sequencing. When quantifying fixed

cells, it is critical to avoid counting cell debris to avoid overestimating the number of cells.


High Quality Sample



Aggregation



Debris

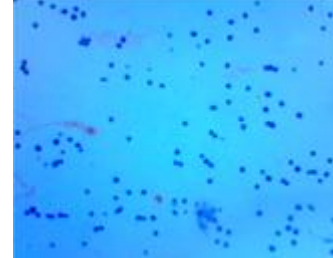

Example trypan blue stained fixed cells

Centrifugation

- A swinging bucket rotor should be used for all high-speed centrifugation steps in this protocol. The use of a fixed-angle rotor will lead to substantial cell loss.

Maximizing Cell Recovery

- It is critical to use the recommended NEB RNase Inhibitor, Murine (New England Biolabs®), in addition to the RNase Inhibitor present in the Cell Fixation kit to maximize RNA integrity and quality.
- It is critical to thoroughly resuspend the cells after centrifugation throughout the protocol. Resuspend by slowly and repeatedly pipetting until no clumps are visible. Ideally this should be verified with microscopy.
- We do not recommend wide bore pipette tips as they make it difficult to resuspend cell pellets adequately.
- Ensure that the 1.5 mL and 15 mL centrifuge tubes are polypropylene, as polystyrene tubes will lead to substantial sample loss.
- When using Evercode Mouse TCR/BCR Fixation kits for the first few times, we recommend retaining the supernatants removed in Section 2 steps 8 and 18 and if relevant, steps 8 and 18 in Appendix B3. In the unlikely event of very high sample loss, these supernatants can be analyzed to identify points for optimization.

Removing Supernatant from Plates

- When using the plate-based Evercode Fixation v3 workflow for the first time, we recommend practicing removing liquid from deep well plates without disturbing the bottom of the well.
- When removing supernatant, keep the pipette tips along the side of the wells to avoid disturbing the pellets. After removing most of the supernatant, we recommend tilting the plate at a 90 degree angle. This makes it easier to see the removal of residual supernatant from bottom and side of the plate.
- Less than 20 μL of supernatant should be left in each well after steps 10 and 26 in Section 2. The figure below shows different volumes of supernatant when viewed from the bottom of the deep well plate.

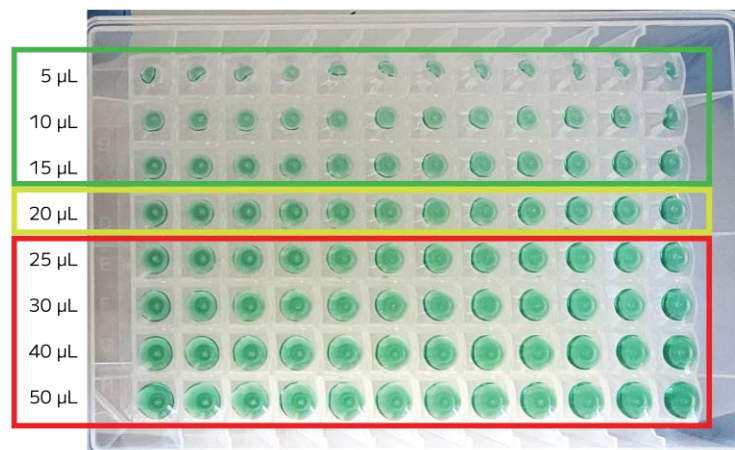


Figure 3: Bottom view of a deep well plate with different volumes of residual supernatant.

Reagent Stability

- Reagents in the Cell Fixation Reagents box should not be frozen and thawed more than 3 times.
- The NEB RNase Inhibitor, Murine, should be stored as recommended by the manufacturer.
- If the kit is going to be used more than 4 times, the reagents should be aliquoted into nuclease-free 1.5 mL tubes and stored at -20°C until use. We do not recommend making single use aliquots to minimize the impact of evaporation during storage.
- To avoid pipetting $<2 \mu\text{L}$ of RNase Inhibitor, we do not recommend preparing less than 2 reactions of Cell Storage Master Mix. In addition, DMSO should be added fresh prior to use in the protocol to a final concentration of 5%. If aliquoting is required, the master

mix can be prepared without DMSO, split into aliquots, and stored at -20°C for up to a month.

Storage of Fixed Samples

- Fixed samples can be stored at -80°C for up to 6 months. Fixed samples should not be refrozen after thawing.
- When possible, we recommend splitting samples into aliquots after fixation in Section 2 Step 32. The aliquots should be at least $20\ \mu\text{L}$ when stored in $0.2\ \text{mL}$ PCR tubes or tube strip(s).
- We recommend making a $20\ \mu\text{L}$ counting aliquot for each sample. This aliquot should be used to update sample concentrations before starting the Evercode workflow. The loss after freeze/thaw should match the sample. This will minimize the time between sample thawing and the start of the Evercode workflow.

Parts List

The Evercode Cell Fixation v3, mid-throughput plate workflow requires Cell Fixation Reagents and Cell Fixation Enhancer boxes.

Cell Fixation Reagents, 12 reactions. Store at -20°C, PN CF100

LABEL	ITEM	PN	FORMAT	QTY
	Prefixation Buffer	CF101	8 mL bottle	1
	Storage Buffer	CF102	2 mL tube	1
	Fixative Solution A	CF103	1.5 mL tube	1
	Fixative Solution B	CF104	1.5 mL tube	1
	Permeabilization Solution	CF105	1.5 mL tube	1
	Fix and Perm Stop Buffer	CF106	8 mL bottle	1
	RNase Inhibitor	CF107	1.5 mL tube	1
	DMSO	CF108	1.5 mL tube	1

Cell Prefixation Enhancer, 12 reactions.

Store at 4°C, PN CF200

LABEL	ITEM	PN	FORMAT	QTY
	Prefixation Enhancer	CF201	1.5 mL tube	1

30 μ M Plate Strainer*. Store at Room Temperature

LABEL	ITEM	PN	FORMAT	QTY
N/A	Adhesive Plate Strainer, 30 μ M	PS030	Plastic sleeve	1

70 μ M Plate Strainer*. Store at Room Temperature

LABEL	ITEM	PN	FORMAT	QTY
N/A	Adhesive Plate Strainer, 70 μ M	PS070	Plastic sleeve	1

100 μ M Plate Strainer*. Store at Room Temperature

LABEL	ITEM	PN	FORMAT	QTY
N/A	Adhesive Plate Strainer, 100 μ M	PS100	Plastic sleeve	1



Note: * Only one mesh size of Plate Strainer is required for the Evercode Cell Fixation v3 workflow. Select an appropriate mesh size for each sample type.

User Supplied Equipment and Consumables

The following materials and equipment are required to perform the protocol but are not provided within the kit. This list does not include standard laboratory equipment, such as freezers.

Reagents

ITEM	SUPPLIER	PN	NOTES
● NEB RNase Inhibitor, Murine	New England Biolabs, Inc.	M0314S M0314L	Choose based on experimental size. Do NOT substitute RNase Inhibitor as it will lead to decrease in sample quality.



CRITICAL! Refer to Section 1.2 of this user guide to ensure you purchase the correct size of ● NEB RNase Inhibitor, Murine from New England Biolabs. For the 48 reactions fixation, it is necessary to use the M0314L reagent.

Equipment

ITEM	SUPPLIER	PN	NOTES
Centrifuge with Swinging Bucket Rotor	Various Suppliers	Varies	Compatible with 96 deep well plates and capable of reaching 4°C.
Microcentrifuge	Various Suppliers	Varies	Compatible with 1.5 mL tubes.
1-channel: P20, P200, P1000 12-channel: P20, P200	Various Suppliers	Varies	Or 8-channel pipettes can be substituted for 12-channel pipettes.
Hemocytometer	Sigma-Aldrich®	Z359629	Or other cell counting device.
Styrofoam Cooler	Various Suppliers	Varies	(Optional) If storing fixed samples before processing with an Evercode Whole Transcriptome kit.
Water bath	Various Suppliers	Varies	(Optional) If preparing aliquots to count the day before running a downstream Evercode kit. Or equivalent thermomixer, heat block, or bead bath capable of holding temperature at 37°C.

Consumables

ITEM	SUPPLIER	PN	NOTES
SealPlate®	Excel Scientific®	100-SEAL-PLT	Or equivalent PCR plate seals. Note that many clear plastic seals are not designed for storage at -80°C, so we recommend using foil plate seals if storing fixed samples in PCR plates.
TempPlate® EXT Sealing Foil	USA Scientific®	2998-0100	(Optional) If storing fixed samples in a PCR plate. Note that many clear plastic seals are not designed for storage at -80°C.
Protein LoBind® Plate	Eppendorf®	951033308 (1 mL) 0030504305 (2 mL)	Or equivalent polypropylene, nuclease-free, v-bottom, 1 mL or 2 mL deep well plates. Do not substitute polystyrene plates as it will lead to substantial cell loss. If possible, we recommend using protein low-binding plates.
Reagent basins	Various Suppliers	Varies	(Optional) If blocking deep well plates with BSA. Sterile, nuclease-free reagent basins.
Eppendorf twin.tec® PCR Plate 96 LoBind®	Eppendorf	0030129504	Or equivalent DNA low-binding, nuclease-free PCR plate capable of holding 270 µL.
Falcon® High Clarity PP Centrifuge Tubes, 15 mL	Corning®	352097	Or equivalent polypropylene centrifuge tubes.
TempAssure® PCR 8-Tube Strips, 0.2 mL	USA Scientific®	1402-4700	Or equivalent nuclease-free 0.2 mL PCR tubes.
RNaseZap™ RNase Decontamination Solution	Thermo Fisher Scientific	AM9780	Or equivalent RNase decontamination solution.
Pipette Tips TR LTS 20 µL, 200 µL, 1,000 µL	Rainin®	17014961, 17014963, 17014967	Or appropriate DNA low-binding, DNase/RNase-free, and filtered pipette tips. Do not use wide bore tips.
Trypan Blue	Various Suppliers	Varies	Or alternative viability dyes, such as AO/PI.

ITEM	SUPPLIER	PN	NOTES
TrypLE™ Express OR TrypLE Select	Thermo Fisher Scientific	12605010 OR 12563011	(Optional) If using adherent cells. Trypsin is not recommended due to variable RNase levels.
Gibco™ Bovine Albumin Fraction V (7.5% solution)	Thermo Fisher Scientific	15260037	(Optional) If blocking tubes with BSA. Chosen due to its low RNase activity. Contact applications support for alternatives.

Section 1: Set Up

1.1. Block Tubes with BSA

Although not required, blocking plates with BSA can increase cell retention. When Protein LoBind plates are not available, we recommend blocking plates, especially for samples with low cell inputs or cells prone to aggregation.

To block tubes:

1. Prepare a fresh 1% BSA as follows, depending on the number of samples being processed.

1% BSA		
Number of Samples	1	12
Nuclease-free water (not supplied)	1.82 mL	21.8 mL
Gibco Bovine Albumin Fraction V (7.5% solution) (not supplied)	280 μ L	3.4 mL
Total Volume	2.1 mL	25.2 mL

2. Fill a new basin with 1% BSA, refilling as needed to complete steps 3 and 4.
3. For each sample, add **1 mL** of 1% BSA to a well of a polypropylene, nuclease-free, v-bottom, 1 mL or 2 mL deep well plate.
4. Repeat step 3 with a second deep well plate.
5. Add new plate seals and invert once to fully coat the wells.
6. Incubate for **30 minutes** at room temperature.
7. Remove the plate seals. Decant and discard the 1% BSA.
8. Add new plate seals.
9. Centrifuge the plates for **1 minute** at 100 x g at room temperature.
10. Remove the plates from the centrifuge and remove the plate seals.
11. Remove any remaining solution from the bottom of the wells with a multichannel P200.

12. Without sealing the plates, air dry for **30 minutes** in a biosafety cabinet at room temperature.
13. Proceed to Section 1.2 or store sealed BSA-coated plates at 4°C for up to 4 weeks.

1.2. Prepare Master Mixes

The reagents listed in this section and the master mix volumes are for the 12-reaction kit. For processing more than 48 samples using the 96-reaction kit, refer to the corresponding reagents and master mix in the High-Throughput Evercode Cell Fixation v3 user manual.

Master mixes are prepared and dispensed into 96 well PCR plates. This makes it easier to transfer reagents into 96 well plates during fixation. Plate Strainers are also pre-cut, if needed. Master mixes should be prepared just prior to fixation.

To prepare master mixes:

1. Fill a bucket with ice. Gather the following items and handle as indicated below.

ITEM	SOURCE	FORMAT	HANDLING AND STORAGE
○ Prefixation Buffer	Cell Fixation Reagents (-20°C)	8 mL bottle	Thaw at room temperature then immediately store on ice. Mix by inverting each tube/bottle. Do not vortex.
● Storage Buffer	Cell Fixation Reagents (-20°C)	2 mL tube	
● Fixative Solution A	Cell Fixation Reagents (-20°C)	1.5 mL tube	
● Fixative Solution B	Cell Fixation Reagents (-20°C)	1.5 mL tube	
● Permeabilization Solution	Cell Fixation Reagents (-20°C)	1.5 mL tube	
○ Fix and Perm Stop Buffer	Cell Fixation Reagents (-20°C)	8 mL bottle	
● DMSO	Cell Fixation Reagents (-20°C)	1.5 mL tube	Thaw and store at room temperature. Mix by inverting the tube.
● RNase Inhibitor	Cell Fixation Reagents (-20°C)	1.5 mL tube	Store on ice immediately before use. Do not vortex.
● NEB RNase Inhibitor, Murine	Purchased Separately	1.5 mL tube	
● Prefixation Enhancer	Cell Prefixation Enhancer (4°C)	1.5 mL tube	

ITEM	SOURCE	FORMAT	HANDLING AND STORAGE
Plate Strainer (30µm, 70µm, or 100µm)	Plate Strainer (Room Temperature)	Plastic sleeve	Keep at room temperature.

2. Prepare the Cell Prefixation Master Mix in a new tube as follows. Mix thoroughly by pipetting and store on ice.

CELL PREFIXATION MASTER MIX			
Number of Samples	1	12	24
○ Prefixation Buffer	203.5 µL	2.44 mL	4.88 mL
● RNase Inhibitor	2.75 µL	33 µL	66 µL
● Prefixation Enhancer	13.75 µL	165 µL	330 µL
Total Volume	220 µL	2.64 mL	5.28 mL



CRITICAL! Reagents in the Cell Fixation Reagents box can be frozen and thawed 3 times. If the kit will be used more than 4 times, aliquots should be made. See Reagent Stability in Important Guidelines for details.

3. For each fixation reaction, dispense **210 µL** of the Cell Prefixation Master Mix into a new 96 well PCR plate. Seal, clearly label and store on ice.



CRITICAL! Do not dispense this master mix into a reagent reservoir because there is not sufficient overage to ensure recovery of the volume required for the fixation. This is true for all the master mixes and reagents dispensed into 96 well PCR plates in this section.

4. Prepare the Cell Fixative Master Mix in a new tube as follows. Mix thoroughly by pipetting and store on ice.

CELL FIXATIVE MASTER MIX			
Number of Samples	1	12	24
● Fixative Solution A	40 µL	480 µL	960 µL
● Fixative Solution B	40 µL	480 µL	960 µL
Total Volume	80 µL	960 µL	1.92 mL

- For each fixation reaction, dispense **77 μL** of the Cell Fixative Master Mix into a new 96 well PCR plate. Seal, clearly label and store on ice.
- Prepare the Cell Storage Master Mix in a new tube as follows. Mix thoroughly by pipetting and store on ice.

CELL STORAGE MASTER MIX			
Number of Samples	1	12	24
● Storage Buffer	112.5 μL	1.35 mL	2.7 mL
● RNase Inhibitor	1.5 μL	18 μL	36 μL
● NEB RNase Inhibitor, Murine	6 μL	72 μL	144 μL
● DMSO	6 μL	72 μL	144 μL
Total Volume	126 μL	1.51 mL	3.024 mL

- For each fixation reaction, dispense **115 μL** of the Cell Storage Master Mix into a new 96 well PCR plate. Seal, clearly label and store on ice.
- For each fixation reaction, dispense **40 μL** of the ● Permeabilization Solution into a new 96 well PCR plate. Seal, clearly label and store on ice.
- Mix the ●OFix and Perm Stop Buffer by inverting the tube 5x. Do not vortex.
- For each fixation reaction, dispense **270 μL** of the ●OFix and Perm Stop Buffer into new 0.2 mL tube strip(s) or 96 well PCR plate. Clearly label and store on ice.
- Each sample is strained twice into 2 different plates in Section 2, so Plate Strainer(s) may need to be cut so they can be applied separately in the two steps with 1 strainer well per sample at both steps. With the Plate Strainers still in the plastic sleeve, cut along the green plastic with sterile scissors, a razor blade, or a scalpel that has been cleaned with RNaseZap.



CRITICAL! Do not directly touch the mesh of Plate Strainer(s) to anything prior to straining.

CRITICAL! After opening their plastic sleeves, Plate Strainers should be stored in their original plastic sleeve until use and used on the same day they were opened.

- Proceed immediately to Section 2.

Section 2: Fixation

2.1. Cell Fixation

After the initial centrifugation to remove the buffer/medium from the single cell suspension, cells are transferred to Cell Prefixation Master Mix. Reagents are added to fix and permeabilize cells, and then stop these reactions. Cells are resuspended in Cell Storage Master Mix and stored at -80°C . This fixation protocol relies on multichannel pipettes and the reagents pre-dispensed 96 well plates in Section 1 to fix up to 48 samples in parallel.

To fix cells:

1. Cool the centrifuge with a swinging bucket rotor to 4°C .
2. Fill a bucket with ice.
3. Prepare a hemocytometer, flow cytometer, or other cell counting device.
4. Count the cells in the single cell suspension with a hemocytometer or alternative counting device and record the count. Keep cells on ice during counting and work quickly to minimize time on ice prior to fixation.
5. With the plate on ice, transfer 100,000 to 1 million cells from each sample into the wells of polypropylene, nuclease-free, v-bottom, 1 mL or 2 mL deep well plate (or BSA-coated deep well plate if prepared in Section 1.1).
6. Add a new plate seal.
7. Centrifuge the plate in a swinging bucket rotor for **5-10 minutes** at 200-500 x g at 4°C .

CRITICAL! Use of a fixed-angle rotor in this protocol will lead to substantial cell loss.



CRITICAL! Ideal centrifugation speed and duration should be determined for each sample type to optimize retention and resuspension efficiencies. See the Important Guidelines section for details.

CRITICAL! Move quickly and handle the samples gently to avoid dislodging the pellet, which will impact data quality.

8. Remove the plate from the centrifuge, remove the plate seal, and store on ice.

9. With a multichannel P200, slowly aspirate and discard all but ~100 μL of supernatant from each well. Keep the pipette tips along the side of the wells to avoid disturbing the pellets.



CRITICAL! Do not reuse any tips across rows throughout this protocol. Never place a tip that has entered one of the wells into a different well or back into a master mix.

10. Tilt the plate 90 degrees. With a multichannel P200, remove and discard the remaining supernatant. Keep the pipette tips along the side of the wells to avoid disturbing the pellets.



CRITICAL! Less than 20 μL of supernatant should be left in each well. See the Important Guidelines section for an example image with varying amounts of residual supernatant.



Note: Tilting the plate 90 degrees makes it easier to visualize removal of residual supernatant.

11. With the plate on ice, fully resuspend each pellet in **180 μL** of Cell Prefixation Master Mix.
12. Apply a Plate Strainer to a new 96 deep well plate (or BSA-coated deep well plate if prepared in Section 1.1) by peeling off the backing, carefully aligning over the wells, and placing on the surface of the plate.



CRITICAL! Do not directly touch the mesh of Plate Strainer(s) to anything prior to straining.



Note: Plate Strainer(s) will need to be cut if processing fewer than 96 samples. See Step 11 in Section 1.2 for details.

13. Without touching the mesh, fully adhere the Plate Strainer by pressing a pipette tip or plate sealer cleaned with RNaseZap along the green plastic on the edges and between the wells of the plate.
14. With both plates on ice, pipette **187.5 μL** of each sample through the strainer into the new 96 deep well plate. Discard the original plate.



Note: To ensure that all of the liquid passes through the strainer, firmly press the tip of the pipette against the filter and steadily depress down the pipette plunger. All of the liquid should pass through the strainer in **~1 second**.

15. Tap the plate 3x on the benchtop to move liquid to the bottom of the wells.
16. Carefully peel off the Plate Strainer and discard.
17. With the plate on ice, add **62.5 μ L** of Cell Fixative Master Mix to each well and mix immediately by pipetting exactly 3x.



CRITICAL! Do not perform additional mixing at this step.

18. Incubate on ice for **10 minutes**.



Note: Start the timer after adding Cell Fixative Master Mix to the first row of the plate.

19. With the plate on ice, add **20 μ L** of ● Permeabilization Solution to each well with a multichannel P20 and mix immediately by pipetting 3x with a multichannel P200 set to 180 μ L.

20. Incubate on ice for **3 minutes**.



Note: Start the timer after adding ● Permeabilization Solution to the first row of the plate. This incubation can be extended by 2 additional minutes up to a total of **5 minutes** without negatively impacting performance.

21. With the plate on ice and changing tips between additions, add **250 μ L** of ○ Fix and Perm Stop Buffer to each well with a multichannel P200 set to 125 μ L and mix immediately by gently pipetting 3x.
22. Add a new plate seal.
23. Centrifuge the plate in a swinging bucket rotor for **5-10 minutes** at 200-500 x g at 4°C.
24. Remove the plate from the centrifuge, remove the plate seal, and store on ice.
25. With a multichannel P200 set to 200 μ L, slowly aspirate and discard **400 μ L** of supernatant from each well. Keep the pipette tips along the side of the wells to avoid disturbing the pellets.

26. Tilt the plate 90 degrees. With a multichannel P200, remove and discard the remaining supernatant from each well. Keep the pipette tips along the side of the wells to avoid disturbing the pellets.



CRITICAL! Less than 20 μL of supernatant should be left in each well after steps 25-26. See the Important Guidelines section for example images of varying volumes of residual supernatant.



Note: Tilting the plate 90 degrees makes it easier to visualize removal of residual supernatant.

27. With the plate on ice, fully resuspend each pellet in **75 μL** of Cell Storage Master Mix.

28. Apply a Plate Strainer to a new 0.2 mL 96 well PCR plate by peeling off the backing, carefully aligning over the wells, and placing on the surface of the plate.



CRITICAL! Do not directly touch the mesh of Plate Strainer(s) to anything prior to straining.

29. Without touching the mesh, fully adhere the Plate Strainer by pressing a pipet tip or plate sealer cleaned with RNaseZap along the green plastic on the edges and between the wells of the plate.

30. With both plates on ice, pipette **85 μL** of each sample through the strainer into the new PCR plate. Discard the original plate.

31. Carefully peel off the Plate Strainer and discard.

32. With a multichannel P20 set to 20 μL , mix the samples by pipetting 5x then immediately transfer **20 μL** to new 0.2 mL tube strip(s). Store on ice.



Note: These aliquots can be thawed and counted separately from the remaining sample for downstream processing with Evercode kits.

33. Seal the original plate with a seal that can withstand storage at -80°C or transfer the remainder of the sample to 0.2 mL tube strip(s).



CRITICAL! Many clear plastic seals are not designed for storage at -80°C , so we recommend using foil plate seals if storing fixed samples in PCR plates.

34. Place the samples in a room temperature styrofoam cooler, close the lid, and store at -80°C to slowly cool the samples.



CRITICAL! Storing samples directly in the freezer without controlled cooling may lead to cell damage and compromise data quality.



Safe stopping point: Samples are stable for up to 6 months at -80°C.

35. The day before running the downstream Evercode kit, thaw the 20 μ L aliquots in a water bath set to 37°C in sets of 2-4. Count the cells in the single cell suspension with a hemocytometer or alternative counting device and record the count. Discard any remaining sample from the counting aliquot.

Appendix: Revision History

Version	Description	Date
1.0	Initial release	Feb 2024
1.1	Updated to accommodate Mouse BCR cell fixation	Dec 2024



parsebiosciences.com

support@parsebiosciences.com

