

User Manual

Version 1.3



# Evercode™ BCR with INTEGRA ASSIST PLUS

## Barcoding Workflow

For use with

ECIB1300 (Human)

ECIB1310 (Mouse)

ECIB1320 (Transgenic  
Mouse)

INTERA ASSIST PLUS

# Legal Notices

This document and its contents are proprietary to Parse Biosciences, Inc. ("Parse Biosciences") and are intended solely for use by its customers in connection with the use of the product(s) described herein and for no other purpose. The products may be used solely FOR RESEARCH PURPOSES, AND MAY NOT BE USED IN ANY DIAGNOSTIC OR THERAPEUTIC USE IN HUMANS OR ANIMALS. This document and its contents shall not be used or distributed for any other purpose and/ or otherwise communicated, disclosed or reproduced in any way whatsoever without the prior written consent of Parse Biosciences.

No rights are granted under this document with respect to any of Parse Biosciences' intellectual property rights. The license to use of any products described herein is subject to a separate written agreement between Parse Biosciences and the applicable user.

The instructions in this document must be strictly and explicitly followed by qualified and properly trained personnel in order to ensure the proper and safe use of the product(s) described herein. Parse Biosciences shall have no liability for any direct, indirect, consequential or incidental damages arising out of any failure to use the product(s) in strict compliance with the terms herein.

This document may contain references to third-party sources of information, hardware or software, products, or services and/ or third-party web sites (collectively "Third Party Information"). Parse Biosciences does not control, and is not responsible for, any Third Party Information. The inclusion of Third Party Information in this document does not imply endorsement by Parse Biosciences of the Third Party Information or the third party in any way.

The product(s) described in this document are provided for one-time use by the purchaser and may not be re-used, refurbished or resold. In addition, such product(s) may not be altered, changed or modified by anyone other than Parse Biosciences and its authorized agents, and Parse Biosciences will not be liable for any such alterations, changes or modifications. The product(s) described herein may be covered by one or more of the following patents:

U.S. Pat. No. 10,900,065

U.S. Pat. No. 11,168,355

U.S. Pat. No. 11,427,856

U.S. Pat. No. 11,634,751

U.S. Pat. No. 11,639,519

U.S. Pat. No. 11,680,283

Patents pending in the U.S. and other countries

PARSE and EVERCODE are trademarks of Parse Biosciences, Inc. The information made available herein is the sole property of Parse Biosciences and its licensors.

Copyright (c) 2024-2025 Parse Biosciences, Inc. All Rights Reserved.

# Table of Contents

Workflow .....	4
Important Guidelines .....	6
Parse Reagents.....	10
Parse-Provided Equipment .....	14
INTEGRA Components.....	14
Equipment .....	15
Consumables .....	16
Section 1: Automation Setup & In Situ Barcoding .....	17
1.1. Sample Normalization.....	17
1.2. Load and Pool Round 1 .....	25
1.3. Round 2 Ligation Preparation .....	29
1.4. Round 2 Ligation .....	32
1.5. Round 3 Ligation Preparation .....	36
1.6. Round 3 Ligation .....	40
1.7. Pre-Lysis.....	43
1.8. Lysis and Sublibrary Generation.....	46
Appendices .....	49
Appendix A: Pipetting Programs .....	49
Appendix B: Troubleshooting.....	57
Appendix C: Revision History.....	58

## Workflow


The Evercode combinatorial barcoding workflow is now compatible with the INTEGRA ASSIST PLUS to enable large-scale single cell RNA-Seq experiments through a robust, semi-automated process.


Evercode split-pool combinatorial barcoding brings a simple workflow to large-scale single cell RNA-seq experiments. The Evercode BCR and Whole Transcriptome kit can profile up to 100,000 cells across up to 48 different biological samples or experimental conditions. Evercode Fixation kits fix and permeabilize cells so they act as individual reaction compartments. This eliminates the need for dedicated microfluidics instrumentation. Through four rounds of barcoding, the transcriptome of each fixed cell is uniquely labeled. The four rounds of barcoding can yield a vast number of possible barcode combinations, which is more than sufficient to uniquely label up to 100,000 cells while avoiding doublets. After sequencing, the Parse Biosciences Analysis Pipeline assigns reads that share the same four barcode combinations to a single cell.

The Evercode BCR and Whole Transcriptome kit is a hybridization-based technology that captures B cells' full length V(D)J sequences.

By integrating the Evercode assay kits with the ASSIST PLUS platform, the semi-automated workflow performs sample normalization and combinatorial barcoding through a series of reverse transcription, ligation, and library preparation steps.

The table below provides a high-level overview of the automated barcoding workflow.

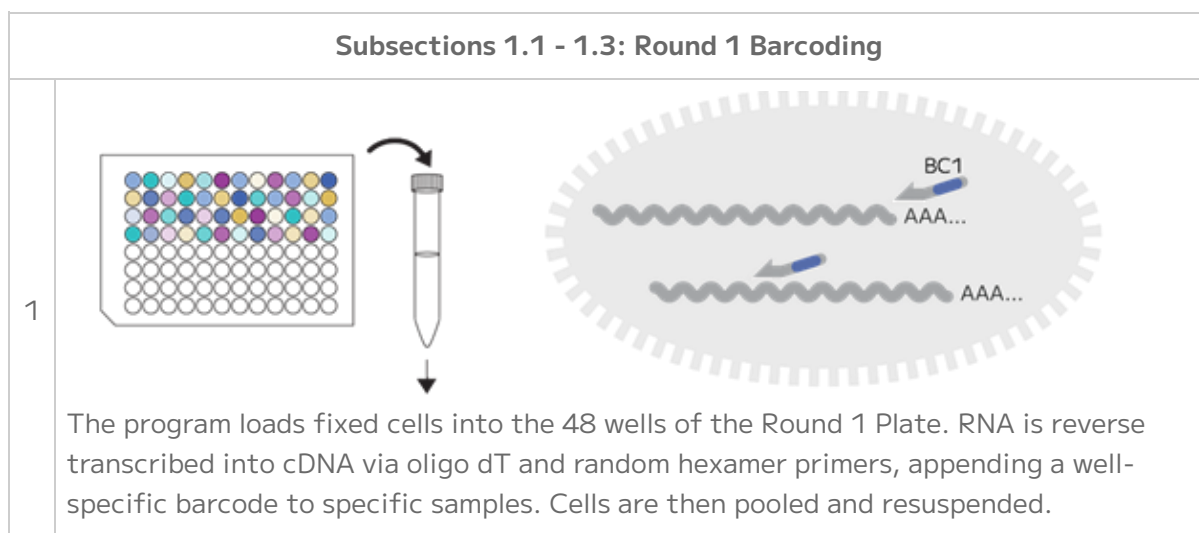

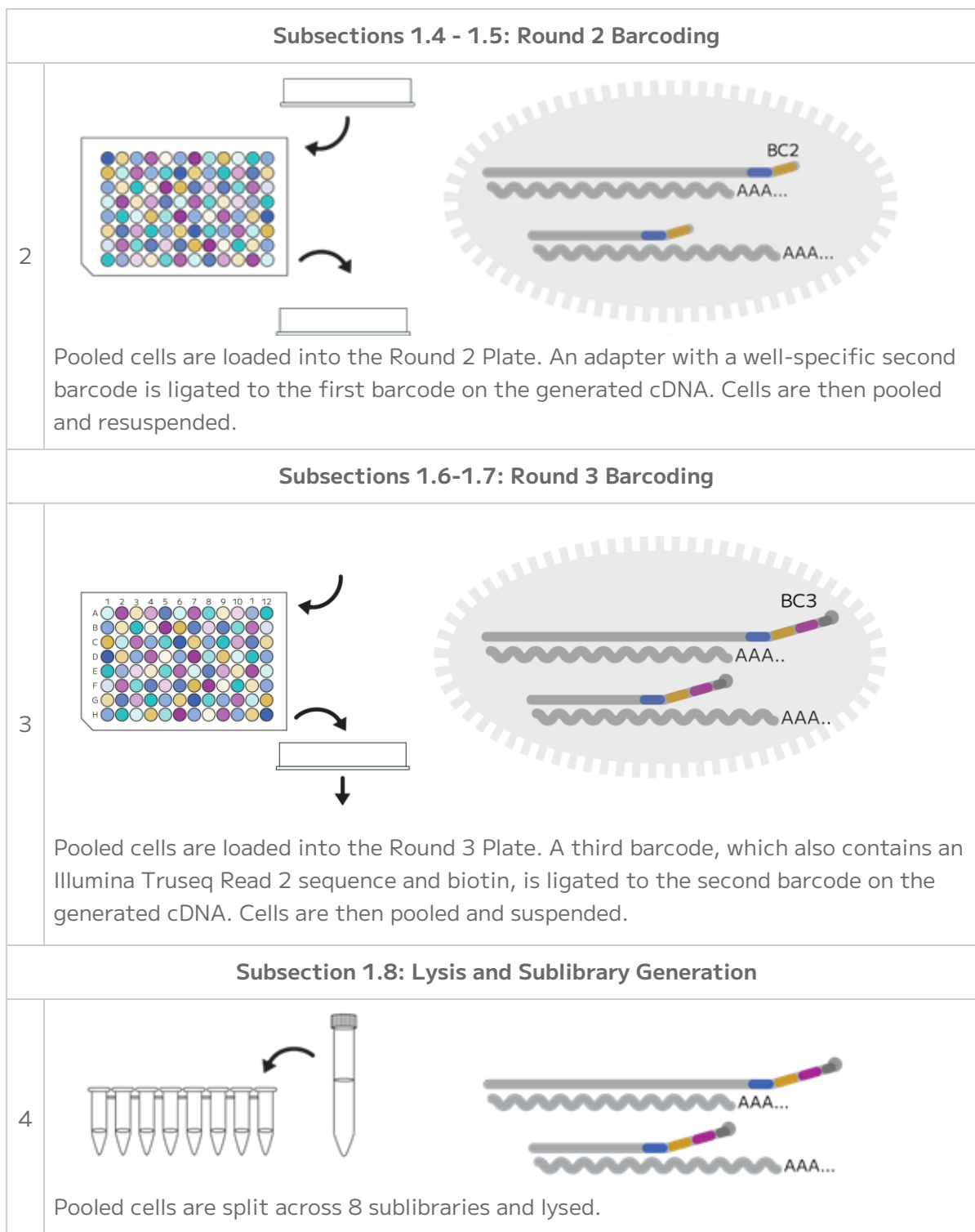

This recommended protocol is intended to be followed when performing the Evercode assays on the INTEGRA ASSIST PLUS liquid handler. The protocol replaces "Section 1. *In Situ* Cell Barcoding", "Section 2. cDNA Purification and Amplification", and "Section 3. Sequencing Library Preparation" of the standard [Evercode BCR User Guide](#).



## Important Guidelines

The following guidelines provide additional information to obtain optimal performance with the Evercode BCR with INTEGRA ASSIST PLUS barcoding workflow.

Comprehensive guidance on optimizing the complete standard Evercode BCR workflow is provided in the Evercode Mega User Guide. For further information on the experimental workflow, please contact [support@parsebiosciences.com](mailto:support@parsebiosciences.com). Please contact [support-us@integra-biosciences.com](mailto:support-us@integra-biosciences.com) for any questions regarding workflow automation or the INTEGRA ASSIST PLUS instrument.

### Sample Input


- This protocol begins with cells previously fixed with an Evercode Cell Fixation v3 kit. It is also compatible with samples fixed with Evercode Cell Fixation v2 kits.
- Even if samples were counted before freezing, we strongly recommend counting cells again after thawing to account for any changes in cell concentration during storage and freeze/thaw. Typically, a 5-15% decrease in concentration after thawing should be expected. These counts are used to determine how cells are loaded in the Round 1 Plate, and counting accuracy at this stage is critical to recover the desired number of cells following barcoding.
- When planning on processing many fixed samples, we recommend setting aside a small "counting aliquot" of each sample at the end of fixation. These counting aliquots can be counted the day prior to starting Section 1 of any Evercode assay kit. The Evercode Fixation User Guides outline recommendations for generating aliquots. Because aliquots have undergone a similar storage time and a freeze/thaw, cell counts from these aliquots will be more representative than using counts from immediately after fixation.
- Aliquots should be thawed in a water bath set to 37°C in sets of 2-4 and counted with a hemocytometer or alternative counting device. Cell counts should be recorded in the Sample Loading Table, and any remaining cell material in the thawed counting aliquot should be discarded.
- Once fixed samples have been thawed, they should not be refrozen.

### Cell Counting and Quality Assessment

- We recommend a hemocytometer for counting, but alternative counting devices can also be used. If possible, validate counts from alternative devices against a hemocytometer when using Evercode kits for the first time.
- When first using Evercode kits, we suggest saving images at each counting step.

- To assess sample quality, we recommend using viability stains like trypan blue or acridine orange and propidium iodide (AO/PI).
- Examples of trypan blue stained fixed cells are shown below. High quality fixed samples have single distinct cells with <5% cell aggregation and no debris. Higher levels of aggregation will lead to elevated doublets after sequencing. When quantifying fixed samples, it is critical to avoid counting cell debris to avoid overestimating the number of cells.


High Quality Sample

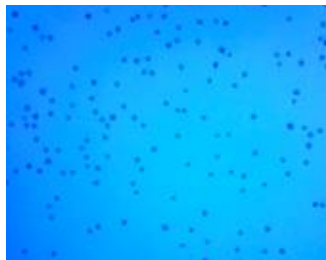

Aggregation



Debris

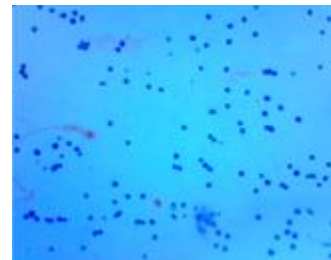

Example trypan blue stained fixed cells.

### Avoiding RNase Contamination

- Take standard precautions to avoid introducing RNases into samples or reagents throughout the workflow. Always wear proper laboratory gloves and use aseptic technique.
- Although RNases are not inactivated by ethanol or isopropanol, they are inactivated by products such as RNaseZap RNase Decontamination Solution (Thermo Fisher Scientific). These can be sprayed on benchtops and pipettes.
- Low retention, filtered pipette tips should be used to reduce RNase contamination from pipettes.

### Centrifugation

- A range of centrifugation speeds and durations are provided in this workflow in order to optimize and balance high sample retention with resuspension efficiency. Comprehensive information to optimize centrifugation conditions for each sample type is provided in any Evercode assay User Guide.
- Use a swinging bucket rotor for all high-speed centrifugation steps in this workflow. A fixed-angle rotor will lead to substantial cell loss.

### Sample Loading Table

- Ensure that the Parse Evercode WT INTEGRA Sample Loading Table being used is the most up-to-date version. The Sample Loading Table can be downloaded from the Parse Biosciences [Customer Support Suite](#). Customer log-in is required to access the Sample Loading Table.
- The Parse Evercode WT INTEGRA Sample Loading Table (Excel spreadsheet) should be completed before starting the experiment.
- If the table isn't working as expected, ensure that Macros are enabled in the Sample Loading Table. Be sure to only edit the colored cells in the table to avoid disturbing any calculations.
- In the unlikely event that samples are not concentrated enough according to the Sample Loading Table (Excel spreadsheet), choose to either:
  - Decrease the "Max number barcoded cells" until the Sample Loading Table no longer gives an error. This will maintain the desired final library proportion between samples but barcode fewer total cells.
  - Proceed as-is. The INTEGRA pipette will dispense the highest possible number of cells based on the existing concentration. This may lead to reduced library yields for samples with particularly low concentrations.

### Plate Sealing

- When sealing or unsealing 96 well plates, do not splash liquid onto the PCR plate seal or between wells. Securing plates in PCR tube racks will minimize this occurrence.
- PCR plate seals may be difficult to remove. Carefully peel off the PCR plate seal while applying downward pressure on the edges of the plate to keep it in the PCR plate rack. Avoid touching the top of the open wells.

### Sample Concentrations

- Dilute the sample with the Pre-Lysis Dilution Buffer to the desired concentration. If the expected sample concentrations are too high, additional Sample Dilution Buffer in the Dilution Accessory Box should be purchased before starting the barcoding workflow.

### PCR Freezer Block

- The PCR freezer block changes color from dark purple to bright pink (too warm) when the block temperature exceeds 7°C. If 70% of the block is bright pink, replace the PCR freezer block with a fully frozen dark purple block.



- To maintain optimal performance and minimize temperature fluctuations, store the PCR Freezer Blocks in a -20°C freezer when not in use.
- Tip pinching may occur when using a fully frozen PCR freezer block for a plate pooling step. To minimize tip pinching, a fully frozen PCR freezer block should be warmed slightly by a few degrees by leaving it at room temperature for 10 minutes before using it on the INTEGRA ASSIST PLUS Deck.
- To ensure a secure positioning within the freezer block, apply firm pressure on the PCR plates until they are securely seated. Avoid touching the top of the open wells.

### INTEGRA ASSIST PLUS Pipetting Programs

- Ensure that the VIALAB program is installed on the same computer that will be used to communicate with the D-ONE Pipetting Module.
- Ensure that the Evercode workflow script precheck has been confirmed prior to running the INTEGRA + Parse workflow. This can be done by running the [Evercode WT with INTEGRA ASSIST PLUS Precheck Scripts](#) available on the Customer Support Suite.
- Ensure that the .csv worklists generated via the Parse Evercode WT INTEGRA Sample Loading Table ("CombinedWTWorksheetYYYYMMDD.csv ") are accessible for upload to the VIALAB program.
- Automated pipetting programs are set up and can be selected using the pipette's user interface. In the INTEGRA + Parse workflow, prompts are built into the programs as checkpoints. When the prompt is followed by a double dash, "--", the program will continue. If the prompt does not contain a double dash, "--", it is followed by another prompt.
- When uploading a new worklist to the VIALAB program, all pipetting settings will be automatically reset to standard default settings and must be readjusted to the correct settings specified for that worklist.












### Deck Loading
















- To prevent malfunctions during the procedure, ensure that all unnecessary labware is cleared from the deck before initiating a new program on the INTEGRA ASSIST PLUS.



## Parse Reagents

The Evercode BCR kit includes the following boxes. Safety Data Sheets for these reagents can be provided upon request.








**-20°C Reagents** Store -20°C, PN HBT100 or MBT100

LABEL	ITEM	PN	FORMAT	QTY
	Round 1 Plate	HBT101 or MBT101	Green semi-skirted 96 well plate	1
	Round 2 Plate	WT102	Blue semi-skirted 96 well plate	1
	Round 3 Plate	WT103	Yellow semi-skirted 96 well plate	1
	Resuspension Buffer	WT104	5 mL tube	1
	Sample Dilution Buffer	WT105	2 mL tube	1
	Round 2 Ligation Buffer	WT106	5 mL tube	1
	Round 2 Ligation Enzyme	WT107	1.5 mL tube	1
	Round 2 Stop Buffer	WT108	2 mL tube	1
	Round 3 Stop Buffer	WT109	5 mL tube	1
	Pre-Lysis Wash Buffer	WT110	5 mL tube	1
	Round 3 Ligation Enzyme	WT111	1.5 mL tube	1



LABEL	ITEM	PN	FORMAT	QTY
	Pre-Lysis Dilution Buffer	WT112	2 mL tube	1
	Lysis Enzyme	WT113	1.5 mL tube	1
	Bead Wash Buffer	WT114	5 mL tube	1
	Wash Buffer 1	WT115	5 mL tube	1
	Wash Buffer 2	WT116	5 mL tube	1
	Capture Enhancer	WT117	1.5 mL tube	1
	Binding Buffer	WT118	1.5 mL tube	1
	Wash Buffer 3	WT119	5 mL tube	1
	Template Switch Buffer	WT120	1.5 mL tube	1
	Template Switch Enzyme	WT121	1.5 mL tube	1
	Template Switch Primer	WT122	1.5 mL tube	1
	cDNA Amp Mix	WT123	1.5 mL tube	1
	cDNA Amp Primers	WT124	1.5 mL tube	1
	Fragm/End Prep Buffer	WT125	1.5.mL tube	1
	Fragm/End Prep Enzymes	WT126	1.5 mL tube	1

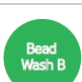LABEL	ITEM	PN	FORMAT	QTY
	Ligation Adapter	WT127	1.5 mL tube	1
	Adapter Ligation Buffer	WT128	1.5 mL tube	1
	Adapter Ligation Enzyme	WT129	1.5 mL tube	1
	Library Amp Mix	WT130	1.5 mL tube	1


**WT + BCR 4°C Reagents.** Store 4°C, PN HBT200 or MBT200

LABEL	ITEM	PN	FORMAT	QTY
	Spin Additive	WT201	1.5 mL tube	1
	Lysis Buffer	WT202	1.5 mL tube	1
	Streptavidin Beads	WT203	1.5 mL tube	1
	Bead Wash Buffer A	GC301	1.5 mL tube	1
	Bead Wash Buffer B	GC302	2 mL tube	1
	Streptavidin Binding Buffer	GC303	5 mL tube	1
	Streptavidin Binder Beads	GC304	1.5 mL tube	1

**BCR -20°C Reagents** Store at -20°C, PN HBT400, MBT400, or MBT500

LABEL	ITEM	PN	FORMAT	QTY
	Hybridization Mix	GC101	1.5 mL tube	1
	Hybridization Enhancer	GC102	1.5 mL tube	1
	Blocker Solution	GC103	1.5 mL tube	1
	Evercode Blocker Solution	GC108	0.5 mL tube	1
 or  or 	Human BCR Panel or Mouse BCR Panel or Transgenic Mouse BCR Panel	GC109 or GC110 or GC111	0.5 mL tube	1
	Enrichment Primer Mix	GC105	1.5 mL tube	1
	Enrichment Amplification Mix	GC106	0.5 mL tube	1

**-20°C Accessory Dilution Buffer** Store at 20°C PN ECAC3901

LABEL	ITEM	PN	FORMAT	QTY
	Sample Dilution Buffer	MG105	2 mL tube	3



**Note:** The Evercode BCR kit ECIB1300 (Human) includes the following boxes: HBT100, HBT200, and HBT400. The Evercode BCR kit ECIB1310 (Mouse) includes the following boxes: MBT100, MBT200, and MBT400.

The Evercode BCR kit ECIB1320 (Transgenic Mouse) includes the following boxes: MBT100, MBT200, and MBT500.

## Parse-Provided Equipment

The following is a list of equipment provided by Parse Biosciences, required to successfully perform the Evercode workflows on the INTEGRA ASSIST PLUS.

ITEM	PN	QTY
Thermochromic PCR Cold Block	NTAC1102	3
Thermochromic PCR Cold Block Riser	NTAC1103	3
Parse Cold Block	NTAC1107	1

## INTEGRA Components

The following are required INTEGRA ASSIST PLUS components needed to run the Evercode assay on the INTEGRA ASSIST PLUS.

ITEM	ITEM TYPE	PN	QTY
Pipette Communication Module for VIAFLO / VOYAGER Pipettes	Accessory	4221	3
ASSIST PLUS Slanted Plate Holder (0°-30°)	Adapter	4510	1
Tip Deck for D-ONE Pipetting Module	Base	4535	1
Dual Reservoir Adapter (ANSI/SLAS footprint)	Adapter	4547	1
ASSIST PLUS Labware Pedestal (+24 mm) for D-ONE, portrait, for 10 mL, 25 mL, 100 mL, reservoirs and plates/reservoirs in SLAS/ANSI format (24 mm Labware Pedestal)	Adapter	4551	1
D-ONE Pipetting Module 1-Ch, 5-1250 µL	Pipette	4532	1
VIAFLO Pipette 12-Ch, 5-125 µL	Pipette	4632	1
VOYAGER Pipette 8-Channel, 5 - 125 µL	Pipette	4722	1



ITEM	ITEM TYPE	PN	QTY
ASSIST PLUS Base Unit	Main	4505	1
Communication/Charging Cable for VIAFLO	Accessory	4226	1
HEATMAG module	Module	4901	1
96 Well Adapter for Magnetic Module	Adapter	4906	1

## Equipment

The following equipment is required to perform the protocol, but are not provided by INTEGRA or by Parse Biosciences. Note that this list does not include standard laboratory equipment, such as freezers.

### Equipment

ITEM	SUPPLIER	PN	NOTES
Microcentrifuge	Various Suppliers	Varies	Compatible with 1.5 mL tubes.
Hemocytometer	Sigma-Aldrich®	Z359629	Or other cell counting device. We recommend validating alternatives relative to a hemocytometer.
Plate Seal Applicator	Various Suppliers	Varies	Capable of adhering plate sealing films to 96 well plates.
T100 Thermal Cycler	Bio-Rad Laboratories®	1861096	Or an equivalent thermocycler compatible with semi-skirted 96 well plates, 0.2mL PCR tubes with up to 100 µL of sample volume, sample heating and cooling from 4-98°C, and a heated lid capable of 30-105°C.
Trypan Blue	Various Suppliers	Varies	Or alternative dyes to assess cell viability, such as AO/PI.

## Consumables

The following is a list of consumables are required to successfully perform Parse Evercode workflows on the INTEGRA ASSIST PLUS.

ITEM	SUPPLIER	PN	QTY
25 mL Basin Reservoir Liner	INTEGRA-Provided	4316	6
Sterilized 40 µm Mini Cell Strainer	DiagnoCine	FNK-HT-AMS-14002	2
10 mL Transport Tube	GlobeScientific	6102S	3
1.5 mL Microtube	Genesee Scientific	21-257	6
PCR Strip Tubes	USA Scientific	1402-4700	8
2 mL Microtubes	Genesee Scientific	21-255	4
8 Row Polystyrene Reservoir	INTEGRA-Provided	6373	2
Semi-skirted 96 well plates	Eppendorf	E951020362	5

## Section 1: Automation Setup & In Situ Barcoding

### 1.1. Sample Normalization

Prior to setting up the Automation workflow, count the cells to assess quality and concentration of the fixed sample(s).

After adjusting the sample(s) to the recommended dilution range, download the Sample Loading Table MACRO (Section 1.1.2), which will be used as reference for allocating the fixed cells into the 96-well PCR plate, ready for the protocol to start.

The program uses the sample dilution buffer on Deck A to normalize fixed samples from plate format on Deck C into intermediate dilution plate on Deck B.

1. Gather the following components and reagents:

ITEM	SOURCE	QTY	HANDLING AND STORAGE
D-ONE Pipetting Module 1-Ch, 5-1250 µL	INTEGRA Component	N/A	
Tip Deck for D-ONE Pipetting Module	INTEGRA Component	N/A	
24 mm Labware Pedestal	INTEGRA Component	N/A	
Thermochromic PCR Cold Block	Parse-Provided	N/A	Pull the Freezer Block with riser from the -20°C freezer and leave them at room temperature for 10 minutes prior to use.
Thermochromic PCR Cold Block Riser	Parse-Provided	N/A	
Parse Cold Block	Parse-Provided	N/A	Keep on ice when not in use.
Semi-Skirted 96 Well PCR Plate	Consumables	N/A	
● Sample Dilution Buffer	-20°C Reagents	1	Thaw at room temperature then store on ice. Mix by inverting 3x.
□ Round 1 Plate	-20°C Reagents	1	Place directly on ice.

2. Download the Parse Biosciences Evercode WT INTEGRA Sample Loading Table. The most current version of the Evercode WT INTEGRA Sample Loading Table can be found on the

Parse Biosciences [Customer Support Suite](#). Customer log-in is required to access the Sample Loading Table.

- Start with the Sample Loading Table tab of the worksheet. Per the instructions in the worksheet, input number of samples (Figure 1).



**Note:** For more information on using the Sample Loading Table, see the "Important Guidelines" section of the User Guide.

This sheet should be filled out prior to starting Section 1.

Step	Instructions
1	Ensure Macros are enabled.
2	Input the number of samples.
3	Input the target number of barcoded cells. <b>Note:</b> The default is 100,000 cells for Evercode WT.
4	Input your sample names.
5	Input the target percentage representation of each sample in the final library. <b>CRITICAL:</b> No percentage can be lower than 2.09%.
If not already done, count the samples as described in Section 1.1 of the Evercode WT User Manual.	
6	Input stock cell concentration for each sample.
7	Prepare the dilutions as described. <b>CRITICAL:</b> Ensure that Sample Dilution Buffer is completely thawed before use.
8	Open the "Plate Configuration" sheet. With the plate on ice, add 14 $\mu$ L of each diluted sample to the appropriate well(s) of the Round 1 Plate as shown in the plate map. <b>CRITICAL:</b> Follow the instructions in the User Guide with respect to sample mixing and changing tips.

**Number of Samples (Step 2):**  **CRITICAL:** We do not recommend editing cells highlighted in grey.

**Target Number Barcoded Cells (Step 3):**

Sample #	Sample Name (Step 4)	Percent of Library (Step 5)	Stock Concentration (cells/ $\mu$ L) (Step 6)	Number of Wells	Targeted Number of Barcoded Cells	Required Sample Concentration (cells/ $\mu$ L)
1		100.00%		48	100000	520
<b>TOTALS:</b>		100.00%		48	100,000	

**Figure 1:** Evercode WT Sample Loading Table.

- While minimizing time on ice, count the number of cells in each sample with a hemocytometer or alternative cell counting device. Record the cell count.
- Enter Target Number of Barcoded Cells, Sample Name, Percent of Library and Stock Concentration for corresponding samples (Figure 2).



**Note:** If any parameter is out of range, the cells will turn pink in color and/or an error message will appear (Figure 2). Be sure to address and adjust worksheet input values appropriately before continuing.

This sheet should be filled out prior to starting Section 1.

Step	Instructions
1	Ensure Macros are enabled.
2	Input the number of samples.
3	Input the target number of barcoded cells. <b>Note:</b> The default is 100,000 cells for Evercode WT.
4	Input your sample names.
5	Input the target percentage representation of each sample in the final library. <b>CRITICAL:</b> No percentage can be lower than 2.09%.
If not already done, count the samples as described in Section 1.1 of the Evercode WT User Manual.	
6	Input stock cell concentration for each sample.
7	Prepare the dilutions as described. <b>CRITICAL:</b> Ensure that Sample Dilution Buffer is completely thawed before use.
8	Open the "Plate Configuration" sheet. With the plate on ice, add 14 $\mu$ L of each diluted sample to the appropriate well(s) of the Round 1 Plate as shown in the plate map. <b>CRITICAL:</b> Follow the instructions in the User Guide with respect to sample mixing and changing tips.

**CRITICAL:** We do not recommend editing cells highlighted in grey.

Number of Samples (Step 2):

Target Number Barcoded Cells (Step 3):

Sample #	Sample Name (Step 4)	Percent of Library (Step 5)	Stock Concentration (cells/ $\mu$ L) (Step 6)	Number of Wells	Targeted Number of Barcoded Cells	Required Sample Concentration (cells/ $\mu$ L)
1	SA-0001	32.00%	500	14	32000	571
2	SA-0002	48.00%	3,300	24	48000	499
3	SA-0003	20.00%	5,250	10	20000	499
<b>TOTALS:</b>		<b>100.00%</b>		<b>48</b>	<b>100,000</b>	


**CRITICAL:** This cell stock concentration is too low.

**Figure 2:** Example error message, noting that the sample stock concentration is too low.

- Navigate to the "INTEGRA Loading Table" tab and check that the Minimum Diluent Needed ( $\mu$ L) does not exceed 1,800  $\mu$ L. If the "Required Number of Sample Dilution Tubes" is greater than 1, additional Sample Dilution Buffer is provided in the Dilution Accessory Box (Figure 3).
- Referring to the "INTEGRA Loading Table" tab, manually load fixed samples into the directed "Sample Location" well into a new Eppendorf semi-skirted plate. Note that any individual sample may require loading into more than one "Sample Location" well. This is the sample stock plate.
- Store the sample stock plate on ice for later use.




**Note:** Example: Sample 1 is loaded into Sample Location A1 and Sample Location A2 with a minimum volume of 131  $\mu$ L. More sample volume (up to a maximum of 200  $\mu$ L total sample volume) can be loaded to reduce bubbles during mixing.

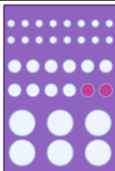Sample Name	Sample Location	Min Sample Stock Needed for Dilution (uL)	*Note: due to semi-skirted plate volumes, multiple wells might be needed for the same samples.	Min Diluent Needed (uL)	Required Number of Sample Dilution Tubes	Sample Dilution Tube Locations
Sample 1	A1	131.0		1955.0	2	
Sample 1	A2	131.0				
Sample 2	A3	84.0				
Sample 3	A4	84.0				
Sample 4	A5	84.0				
	A6					
	A7					
	A8					
	A9					
	A10					
	A11					
	A12					
	B1					

**Figure 3:** INTEGRA Loading Table tab noting Minimum Diluent Needed and Sample Locations.

- Navigate to the "Diluent Volumes" tab and click on the purple box labeled "Generate a Worklist for Import into VIALAB". Save the generated CSV file (called "CombinedWTWorksheet.csv") for later use (Figure 4).

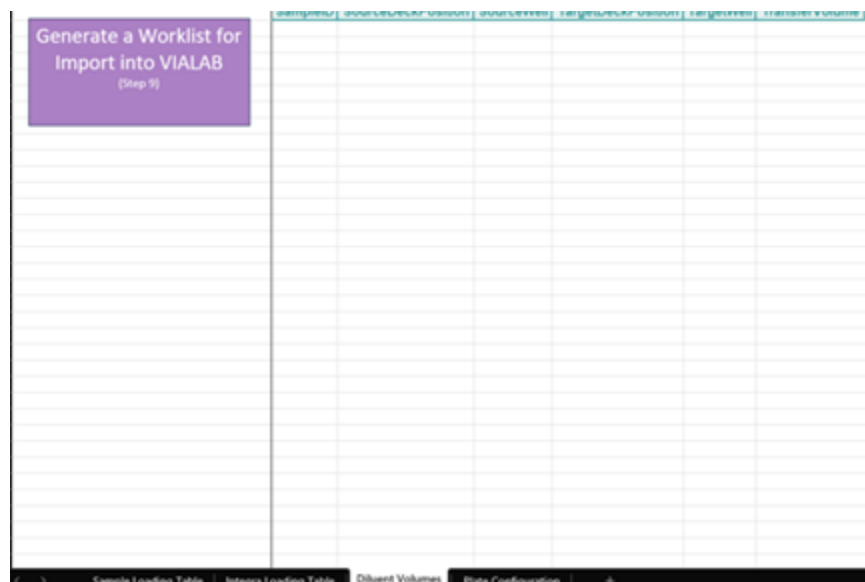

**Figure 4:** Diluent Buffer Volumes tab for generating a VIALAB worklist.

- Navigate to the "Plate Configuration" tab to visualize the final sample location and orientation within the 96 well plate format (Figure 5).

Parse Biosciences  
700 Dexter Ave  
Suite 600  
Seattle, WA 98109



### Evercode WT Round 1 Barcoding Plate Configuration

- For more details on using this Sample Loading Table, see the Important Guidelines section of the User Manual

Use the following plate layout to load samples into the Round 1 Plate in Section 1.2.


	1	2	3	4	5	6	7	8	9	10	11	12
A	1	1	1	1	1	1	1	1	1	1	1	1
B	1	1	2	2	2	2	2	2	2	2	2	2
C	2	2	2	2	2	2	2	2	2	2	2	2
D	2	2	3	3	3	3	3	3	3	3	3	3
E	Do not add samples to these wells.											
F												
G												
H												

Sample Number	Sample Name	Percent Contributing
1	SA-0001	32.00%
2	SA-0002	48.00%
3	SA-0003	20.00%



**Figure 5:** Plate Configuration tab visualizes the sample locations and orientations.


11. Open the VIALAB program **WT S1 St1 DONE V3** and navigate to the “Method” section. In the “O2 Worklist”, under the “Worklist and Volumes” tab, upload the “CombinedWTWorksheet.csv” worklist file generated in Step 9 using the “Import” button (Figure 6).

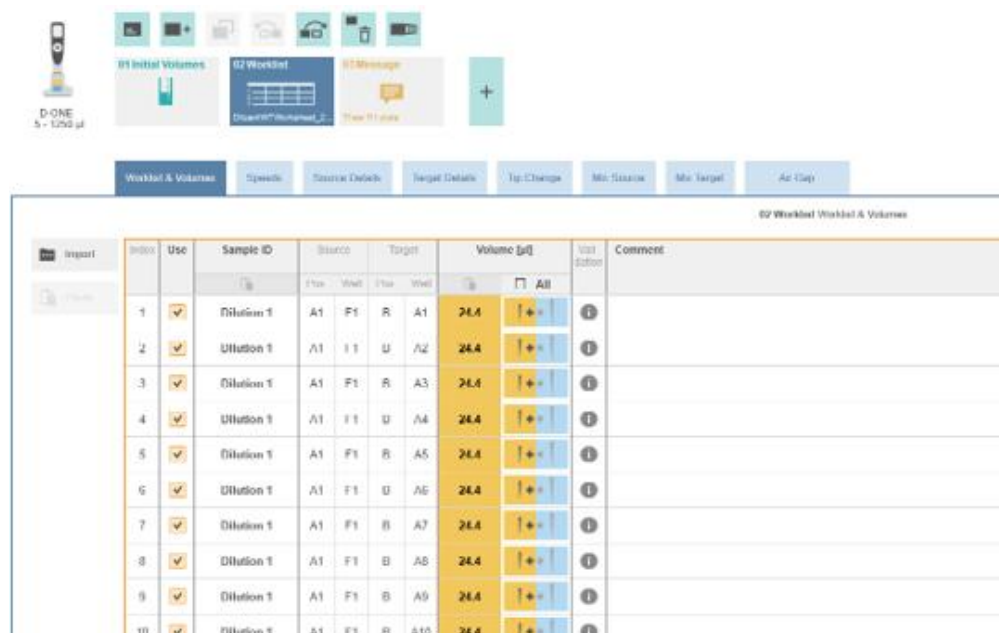

Index	Use	Sample ID	Source	Target	Volume [µL]	Comment
1	✓	Dilution 1	A1 F1 R	A1	24.4	
2	✓	Dilution 1	A1 F1 U	A2	24.4	
3	✓	Dilution 1	A1 F1 R	A3	24.4	
4	✓	Dilution 1	A1 F1 U	A4	24.4	
5	✓	Dilution 1	A1 F1 R	A5	24.4	
6	✓	Dilution 1	A1 F1 U	A6	24.4	
7	✓	Dilution 1	A1 F1 R	A7	24.4	
8	✓	Dilution 1	A1 F1 U	A8	24.4	
9	✓	Dilution 1	A1 F1 R	A9	24.4	
10	✓	Dilution 1	A1 F1 U	A10	24.4	


**Figure 6:** VIALAB worklist generation for diluent volumes using O2 Worklist.

12. Save and transfer your program to your D-ONE Pipetting Module (1-Ch, 5-1250 µL) as follows:
  - a. Connect your computer with D-ONE Pipette using Communication/Charging Cable for VIAFLO.
  - b. Follow the instructions on the pipette menu, as shown in the diagram below:

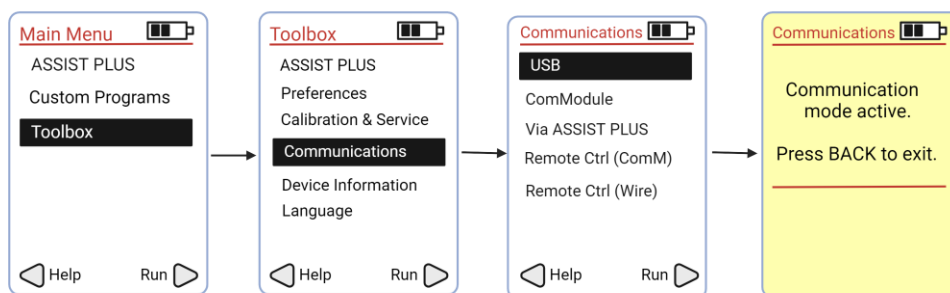

- c. In the VIALAB on your computer, select “Transfer”.

- d. Confirm that the pipette has been recognized by the software. Click "Send to pipette". This will upload the **WT S1 St1 V3** program to the D-ONE Pipette
- e. Disconnect the cable and proceed with the protocol on the INTEGRA ASSIST PLUS platform.
- f. If done correctly, a program named **WT S1 St1 V3** will be found on your pipette as shown in the diagram below.

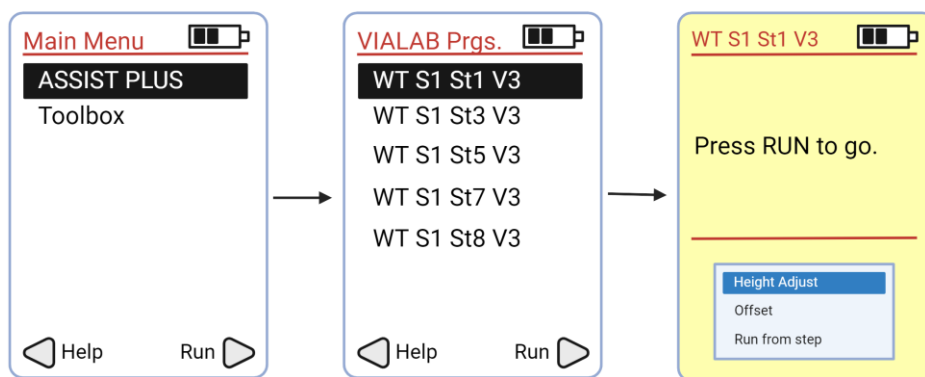

13. Place the Barcoding Reagents in an ice bucket.
14. Remove both Thermochromic PCR Cold Blocks with Thermochromic PCR Cold Block Risers from -20°C and thaw at room temperature for **10 minutes**. Ensure that the corner sizes align between the colored Thermochromic PCR Cold Block and the gray Thermochromic PCR Cold Block Riser (Figure 7).



**Figure 7:** Corners are aligned between the Thermochromic PCR Cold Block and the Thermochromic PCR Cold Block Riser.


15. Set up the 3 Position Universal Deck according to the deck configuration below.



**Note:** Refer to the 'Integra Loading Table' tab of the Parse Biosciences Evercode WT INTEGRA Sample Loading Table. Place the dilution tubes according to the locations

highlighted in pink under 'Sample Dilution Tube Locations' on the deck configuration (Figure 3)

## Deck Configuration

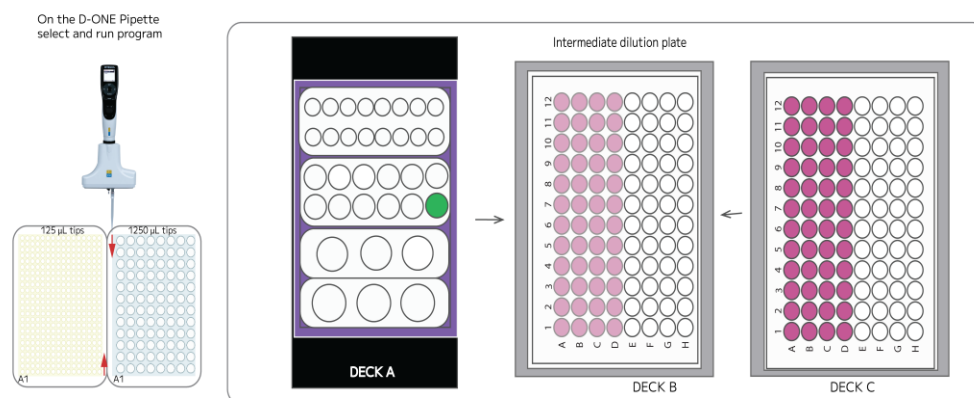

	DECK A	DECK B	DECK C
HARDWARE	Parse Cold Block 24 mm Labware Pedestal	Thermochromic PCR Cold Block Thermochromic PCR Cold Block Riser	Thermochromic PCR Cold Block Thermochromic PCR Cold Block Riser
CONSUMABLES		Semi skirted plate (empty)	Semi skirted plate (with sample)
REAGENTS	● Sample Dilution Buffer		● Samples

16. Attach D-ONE Pipetting Module 1-Ch, 5-1250 µL and the corresponding Tip Deck to the ASSIST Plus.



**Note:** Ensure the Tip Deck and Tip Boxes are securely positioned by pressing them down firmly until they click into place.

17. Remove the reagent caps, then select and run the program **WT S1 St1 V3**.



18. **When prompted**, thaw the  Round 1 Plate using the thermocycling program below.

THAW ROUND 1 PLATE		
Run Time	Lid Temperature	Sample Volume
10 min	70°C	26 µL
Step	Time	Temperature
1	10 min	25°C
2	Hold	4°C

19. Remove Thermochromic PCR Cold Block from -20°C freezer and thaw it at room temperature for the duration of the thermocycling program.

20. Gently remove the  Round 1 Plate from the thermocycler and centrifuge for **1 minute** at 100 x g at 4°C.

## 1.2. Load and Pool Round 1

The program loads the normalized cells from Section 1.1 on Deck B into the  Round 1 Plate on Deck C. After Barcoding Round 1 incubation, move the  Round 1 plate onto Deck B.

The program then pools all the samples in the  Round 1 Plate into row A.

### To load the sample(s):

1. Gather the following components and reagents:



ITEM	SOURCE	QTY	HANDLING AND STORAGE
VIAFLO Pipette 12-Ch, 5-125 µL	INTEGRA Component	N/A	
Thermochromic PCR Cold Block	Parse-Provided	N/A	Pull the Freezer Block with riser from the -20°C freezer and leave them at room temperature for 10 minutes prior to use.
Thermochromic PCR Cold Block Riser	Parse-Provided	N/A	

2. Remove D-ONE Pipetting Module 1-Ch, 5-1250 µL and corresponding Tip Deck. Replace it with VIAFLO 12-Ch 5-125 µL Pipette and corresponding Tip Deck.




**Note:** Before removing the D-ONE Pipetting Module 1-Ch, 5-1250 µL, ensure it is disconnected by pressing the back button until the main menu is displayed.

**Note:** Ensure the Tip Deck and Tip Boxes are securely positioned by pressing them down firmly until they click into place. Remove the tip box lid prior to starting the program.


3. Remove all items from Deck A.
4. Remove the  Round 1 Plate from the centrifuge, place on a stable surface and remove the plate seal.
5. Replace the Thermochromic PCR Cold Block on Deck C with the new Thermochromic PCR Cold Block thawed in Section 1.1.20. Place this old Thermochromic PCR Cold Block back in the -20°C freezer for later use.
6. **When prompted**, place the  Round 1 Plate in the Thermochromic PCR Cold Block located on Deck C. Ensure A1 is oriented towards the bottom left corner. Deck should correspond to the diagram below.

## Deck Configuration

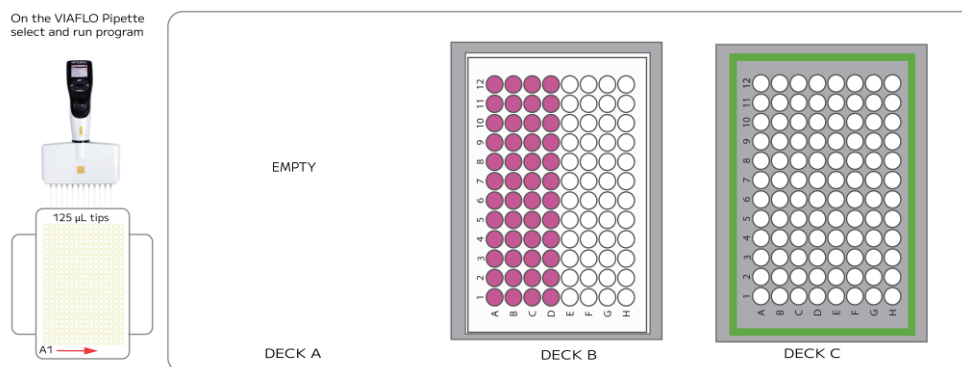

	DECK A	DECK B	DECK C
HARDWARE		Thermochromic PCR Cold Block Thermochromic PCR Cold Block Riser	Thermochromic PCR Cold Block Thermochromic PCR Cold Block Riser
CONSUMABLES			
REAGENTS		● Samples	□ Round 1 Plate

- On the VIAFLO Pipette 12-Ch, select and run the program **WT S1 St2 V3** following the diagram below.

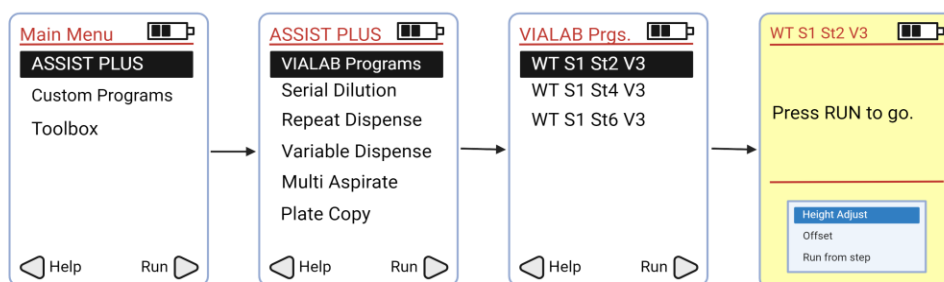

- When prompted**, seal the □ Round 1 Plate from Deck C using a new plate seal. This is best achieved while the plate is secured in a PCR plate rack and on a flat surface.




9. Place the  Round 1 Plate into a thermocycler and run the following program.


BARCODING ROUND 1			
Run Time	Lid Temperature	Sample Volume	
40 min	70°C	40 µL	
Step	Time	Temperature	Cycles
1	10 min	50°C	1
2	12 s	8°C	3
3	45 s	15°C	
4	45 s	20°C	
5	30 s	30°C	
6	2 min	42°C	
7	3 min	50°C	
8	5 min	50°C	1
9	Hold	4°C	Hold

10. Discard the used semi-skirted plate on Deck B. Freeze the Thermochromic PCR Cold Block on Deck B in a -20°C freezer.

11. Move the Thermochromic PCR Cold Block from Deck C to Deck B.

12. **When prompted**, once Barcoding Round 1 thermocycling program is over, place  Round 1 Plate in Thermochromic PCR Cold Block located on Deck B with A1 oriented towards the bottom left corner. Deck layout should correspond to the configuration below.

## Deck Configuration

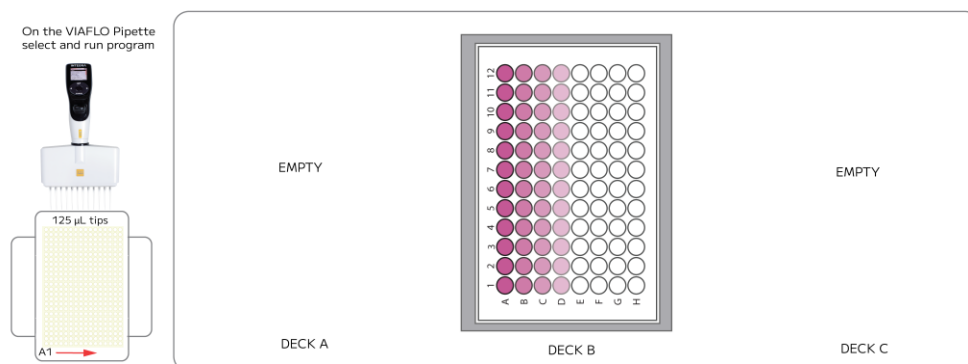




13. Remove the plate seal.

14. Press "Run" on the pipette to continue the program.

15. At the conclusion of the run, we recommend removing the Thermochromic PCR Cold Block with the Thermochromic PCR Cold Block Riser from the deck and storing them in the freezer.

### 1.3. Round 2 Ligation Preparation

The program pools all the samples from the plate on Deck B into a 10 mL tube on Deck C. After centrifugation, removes supernatant and resuspends the cells in  Resuspension Buffer. Mixes  Round 2 Ligation Enzyme and  Round 2 Ligation Buffer with the resuspended cells. The cells are then put on the left basin (A1) on Deck A.

1. Gather the following components and reagents:

ITEM	SOURCE	QTY	HANDLING AND STORAGE
D-ONE Pipetting Module 1-Ch, 5-1250 µL	INTEGRA Components	N/A	
Thermochromic PCR Cold Block	Parse-Provided	N/A	Pull the freezer block with riser from the -20°C freezer and leave them at room temperature for 10 minutes prior to use.
Thermochromic PCR Cold Block Riser	Parse-Provided	N/A	
Parse Cold Block	Parse-Provided	N/A	Keep on ice when not in use.
Tip Deck for D-ONE Pipetting Module	INTEGRA Components	N/A	
10 mL transport tube	Consumables	1	
Dual Reservoir Adapter	INTEGRA Components	N/A	
25 mL Basin Reservoir Liners	INTEGRA-Provided	2	
 Round 2 Plate	-20°C Reagents	1	Place directly on ice.
 Round 2 Ligation Enzyme	-20°C Reagents	1	Place directly on ice. Briefly centrifuge before use.
 Round 2 Ligation Buffer	-20°C Reagents	1	Thaw at room temperature then store on ice. Mix by inverting 3x.
 Resuspension Buffer	-20°C Reagents	1	
 Spin Additive	4°C Reagents	1	Briefly centrifuge. Keep at room temperature.


2. Place the Parse Cold Block on Deck C.

- Place the Dual Reservoir Adapter (INTEGRA logo oriented to the front) on Deck A lined with two new 25 mL basin reservoir liners.
- Configure the deck layout as follows:



**Note:** Centrifuge reagents before loading on Deck C.

### Deck Configuration

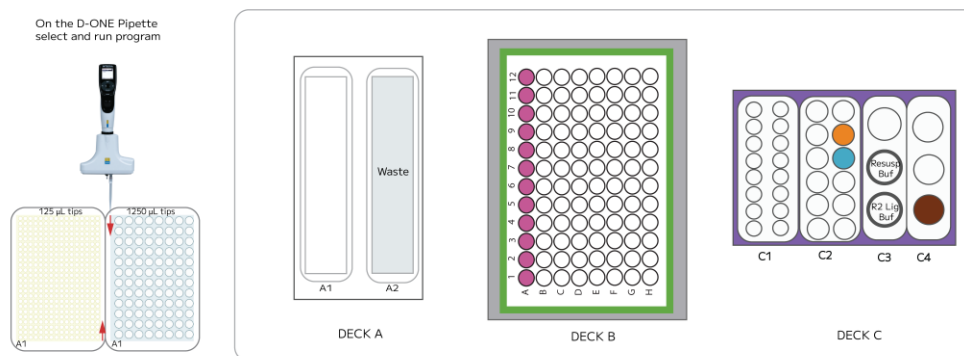

	DECK A	DECK B	DECK C
HARDWARE	Dual Reservoir Adapter	Thermochromic PCR Cold Block Thermochromic PCR Cold Block Riser	Parse Cold Block
CONSUMABLES	25 mL basin reservoir liners		● 10 mL transport tube
REAGENTS	● Samples	□ Round 1 Plate	● Spin Additive ● Round 2 Ligation Enzyme ○ Resuspension Buffer ○ Round 2 Ligation Buffer


- Remove VIAFLO Pipette 12-Ch, 5-125 µL and corresponding Tip Deck. Replace it with the D-ONE Pipetting Module 1-Ch, 5-1250 µL and corresponding Tip Deck.



**Note:** Before removing the VIAFLO Pipette 12-Ch, 5-125 µL, ensure it is disconnected by pressing the back button until the main menu is displayed.

**Note:** Ensure the Tip Deck and Tip Boxes are securely positioned by pressing them down firmly until they click into place. Remove the tip box lid prior to starting the program.

- Remove the reagent caps, then run the program **WT S1 St3 V3** following the diagram below.

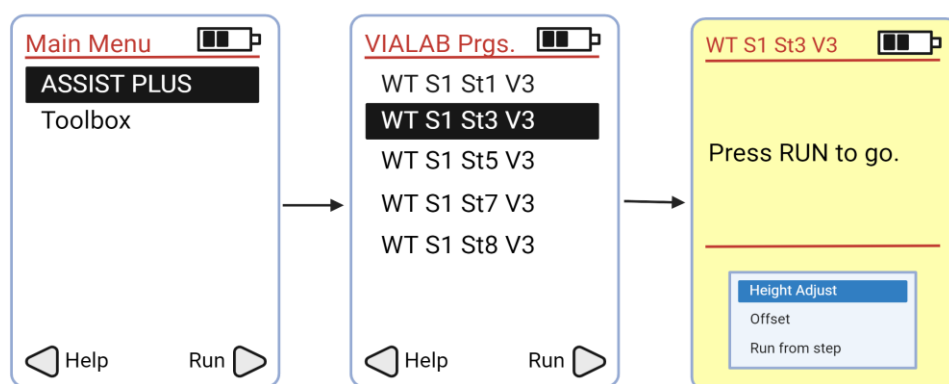

- When prompted**, cap and invert the 10 mL transport tube containing the pooled cells. Centrifuge the sample in a swinging bucket centrifuge cooled to 4°C for **10 minutes** at 250 x g.
- During the 10 minutes spin, **when prompted**, thaw the ☐ Round 2 Plate using the program below for later use. Proceed to the next step while the program is still running.


THAW ROUND 2 PLATE		
Run Time	Lid Temperature	Sample Volume
10 min	70°C	10 µL
Step	Time	Temperature
1	10 min	25°C
2	Hold	4°C

- Remove a new Thermochromic PCR Cold Block with Thermochromic PCR Cold Block Riser from -20°C freezer and thaw for **10 minutes** during centrifugation for later use.
- Once centrifugation is complete, **when prompted**, remove the cap from the 10 mL transport tube and place it back in its original position within the Parse Cold Block 1 located on Deck C4. Immediately proceed to the next step.


11. Press “Run” to continue.

12. Clear Decks B and C. Discard the right basin lines and Deck A. Place the Parse Cold Block on ice.

## 1.4. Round 2 Ligation

The program transfers Cell Suspension Mix from the left reservoir (A1) on Deck A to  Round 2 Plate on Deck B. Following Round 2 Incubation, the program will load Round 2 Stop Buffer in the right reservoir (A2) into each well. After Round 2 Stop incubation, it pools all samples together in the left reservoir (A1) on Deck A.

1. Gather the following items and handle as indicated below:

ITEM	SOURCE	QTY	HANDLING AND STORAGE
VIAFLO Pipette 12-Ch, 5-125 µL	INTEGRA Components	N/A	
Tip Deck for VIAFLO Pipette 12-Ch, 5-125 µL	INTEGRA Components	N/A	
Thermochromic PCR Cold Block	Parse-Provided	N/A	Pull the Freezer Block with riser from the -20°C freezer and leave them at room temperature for 10 minutes prior to use.
Thermochromic PCR Cold Block Riser	Parse-Provided	N/A	
INTEGRA Dual 25mL Basin Reservoir Adapter	INTEGRA-Provided	N/A	
25 mL Basin Reservoir Liners	INTEGRA-Provided	2	
 Round 2 Stop Buffer	-20°C Reagents (Parse Reagents)	1	Thaw at room temperature then store on ice. Mix by vortexing.




2. Remove D-ONE Pipetting Module 1-Ch, 5-1250 µL and corresponding Tip Deck. Replace it with VIAFLO 12-Ch 5-125 µL and corresponding 125 µL Tip Deck.




**Note:** Before removing the D-ONE Pipetting Module 1-Ch, 5-1250 µL, ensure it is disconnected by pressing the back button until the main menu is displayed.




**Note:** Ensure the Tip Deck and Tip Boxes are securely positioned by pressing them down firmly until they click into place. Remove the tip box lid prior to starting the program.


3. Remove the  Round 2 Plate from the thermocycler and centrifuge for **1 minute** at 100 x g at 4°C.
4. While secured in a PCR plate rack on a flat surface, remove the plate seal from the  Round 2 Plate and store on ice.
5. Place the Thermochromic PCR Cold Block thawed during Section 1.3.9 on Deck B.
6. Place the  Round 2 Plate on Deck B in the Thermochromic PCR Cold Block with Thermochromic PCR Cold Block Riser with A1 oriented towards the bottom left corner and remove the seal. The deck should correspond to the configuration below.

## Deck Configuration

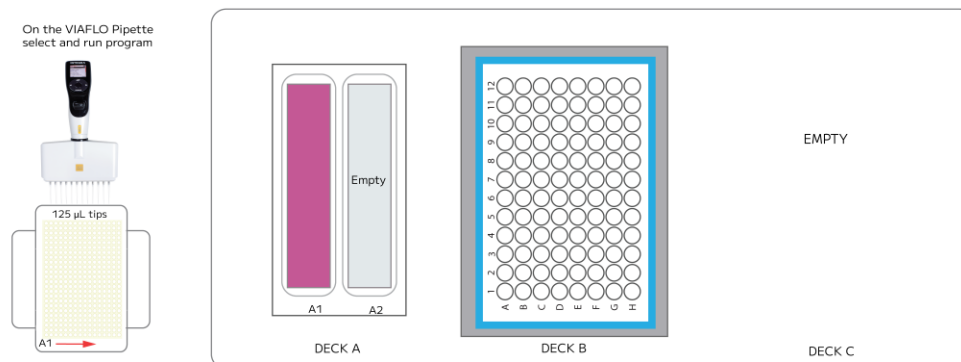

	DECK A	DECK B	DECK C
HARDWARE	Dual Reservoir Adapter	Thermochromic PCR Cold Block Thermochromic PCR Cold Block Riser	
CONSUMABLES	25 mL Reservoir liners		
REAGENTS	● Samples	 Round 2 Plate	

7. Select and run the program **WT S1 St4 V3** following the diagram below.

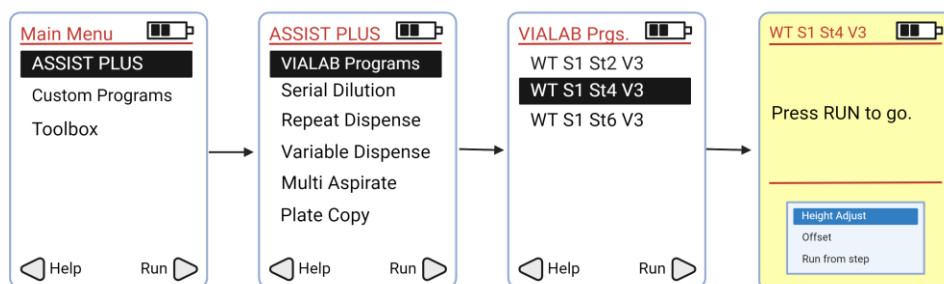


8. **When prompted**, reseal the  Round 2 Plate with an adhesive seal, place it into a thermocycler, and run the following program. Upon completion, proceed immediately to the next step.

BARCODING ROUND 2		
Run Time	Lid Temperature	Sample Volume
15 min	50°C	50 µL
Step	Time	Temperature
1	15 min	16°C
2	Hold	4°C

9. Remove the  Round 2 Plate from the thermocycler and place in a PCR plate rack.
10. Remove the plate seal and place the  Round 2 Plate back on Deck B with A1 oriented towards the lower left corner.
11. **When prompted**, replace the 25 mL basin reservoir liner on the right with a new 25 mL basin reservoir liner.
12. Briefly vortex (2-3 seconds) and centrifuge the Parse ● Round 2 Stop Buffer. **When prompted**, using a P1000 pipette, add the total volume (~1.4 mL) to the right basin (A2) on Deck A.

## Deck Configuration

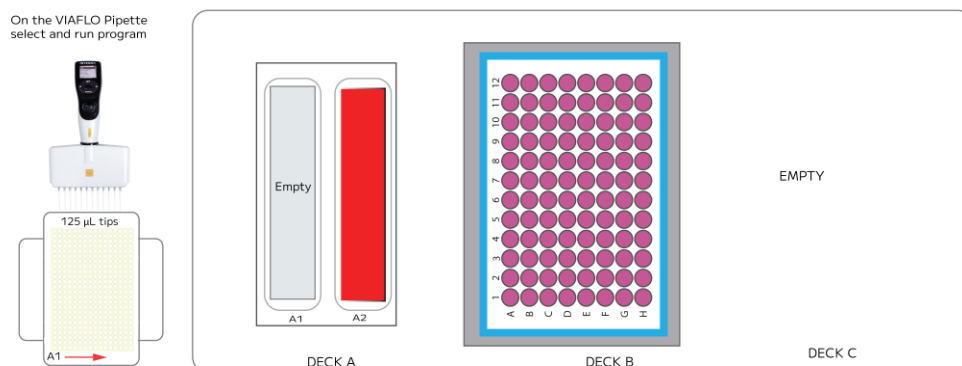


	DECK A	DECK B	DECK C
HARDWARE	Dual Reservoir Adapter	Thermochromic PCR Cold Block Thermochromic PCR Cold Block Riser	
CONSUMABLES	25 mL reservoir liners		
REAGENTS	● Round 2 Stop Buffer	□ Round 2 Plate	




**CRITICAL!** When adding the ● Round 2 Stop Buffer to the reservoir ensure the volume is evenly distributed for optimal pipetting.

- Remove the reagent caps, then press "Run" to continue.
- When prompted**, reseal the □ Round 2 Plate with an adhesive seal, and incubate the □ Round 2 Plate in a thermocycler using the following protocol. Upon completion, proceed immediately to the next step.

ROUND 2 STOP		
Run Time	Lid Temperature	Sample Volume
5 min	50°C	60 µL
Step	Time	Temperature
1	5 min	16°C
2	Hold	4°C



15. Place the  Round 2 Plate within the Thermochromic PCR Cold Block with Thermochromic PCR Cold Block Riser located on Deck B with A1 oriented towards the bottom left.
16. Remove the seal and press "Run".
17. **When prompted**, replace the right basin liner (A2) on Deck A.
18. At the conclusion of the run, we recommend removing the Thermochromic PCR Cold Block with the Thermochromic PCR Cold Block Riser from the deck and storing them in the freezer.

## 1.5. Round 3 Ligation Preparation

The pooled cell suspension in the left reservoir (A1) on Deck A is strained into the 10 mL transport tube on Deck C4. The program adds  Round 3 Ligation Enzyme into the cells and mixes. It then transfers the cell suspension mix to the top reservoir (A1) within the slanted plate holder on Deck B.

1. Gather the following components and reagents:


ITEM	SOURCE	QTY	HANDLING AND STORAGE
D-ONE Pipetting Module 1-Ch, 5-1250 µL	INTEGRA Component	N/A	
Tip Deck for D-ONE Pipetting Module	INTEGRA Component	N/A	
Parse Cold Block	Parse-Provided	N/A	Keep on ice when not in use.
Dual Reservoir Adapter (ANSI/SLAS footprint)	INTEGRA Component	N/A	
ASSIST PLUS Slanted Plate Holder	INTEGRA Component	N/A	
10 mL transport tube	Consumables	1	
25 mL Basin Reservoir Liners	INTEGRA-Provided	2	
40 µm cell strainer	Consumables	1	
125 µL Tip Rack	INTEGRA-Provided	1	
1250 µL Tip Rack	INTEGRA-Provided	1	

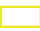ITEM	SOURCE	QTY	HANDLING AND STORAGE
 Round 3 Plate	-20°C Reagents	1	Place directly on ice.
 Round 3 Ligation Enzyme	-20°C Reagents	1	Place directly on ice. Briefly centrifuge before use.

2. Add the ASSIST PLUS Slanted Plate Holder on Deck B, incorporating a **10° tilt**.



**Note:** Ensure that the lowest side is positioned adjacent to Deck A.


3. Place the Parse Cold Block 1 on Deck C.
4. Briefly centrifuge and insert the  Round 3 Ligation Enzyme tube in the appropriate location in the Parse Cold Block.
5. Place a 10 mL transport tube with the 40 µm cell strainer within its respective location in the Parse Cold Block 1 located in the Reagent Block in C4 position. Deck layout should correspond to the Deck Configuration below.

## Deck Configuration




**Note:** When moving the dual reservoir adapter onto the slanted plate holder on Deck B, ensure the INTEGRA logo is oriented towards the front left.

**Note:** Use extra care when moving the cell suspension to avoid spills.

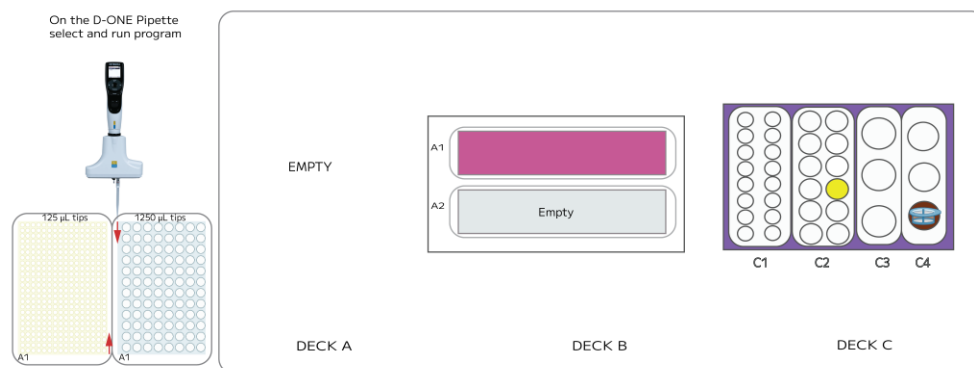

### Slanted Plate Holder (10°) front view

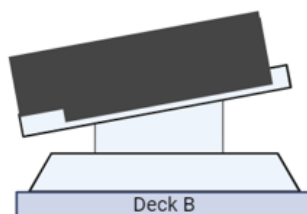

	DECK A	DECK B	DECK C
HARDWARE		Slanted Plate Holder (10°) Dual Reservoir Adapter	Parse Cold Block
CONSUMABLES		25 mL basin reservoir liners	● 10 mL transport tube with cell strainer
REAGENTS		● Samples	● Round 3 Ligation Enzyme


- Remove VIAFLO Pipette 12-Ch, 5-125  $\mu$ L and corresponding tip deck. Replace it with D-ONE Pipetting Module 1-Ch, 5-1250  $\mu$ L and corresponding Tip Deck.




**Note:** Before removing the VIAFLO Pipette 12-Ch, 5-125  $\mu$ L, ensure it is disconnected by pressing the back button until the main menu is displayed.

**Note:** Ensure the Tip Deck and Tip Boxes are securely positioned by pressing them down firmly until they click into place.

- Remove the reagent caps, then select and run the program **WT S1 St5 V3** following the diagram below.

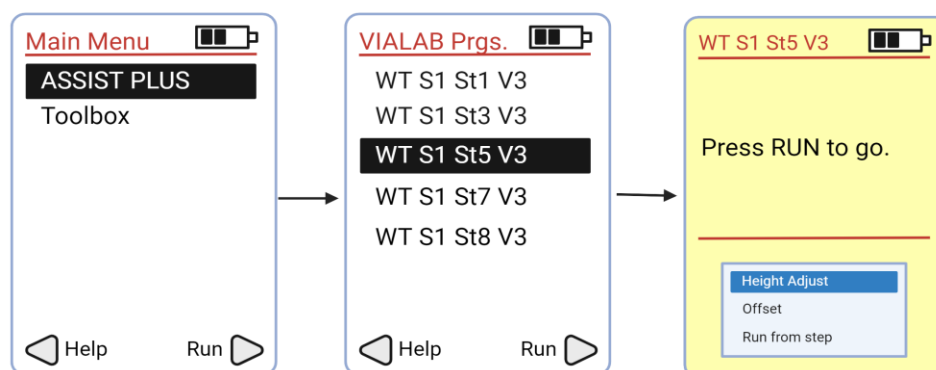

- Proceed to the next step while the program is still running.
- When prompted**, place the  Round 3 Plate into a second thermocycler and run the program below. Proceed to the next step while the program is still running.

THAW ROUND 3 PLATE		
Run Time	Lid Temperature	Sample Volume
10 min	70°C	10 µL
Step	Time	Temperature
1	10 min	25°C
2	Hold	4°C

- Take a new Thermochromic PCR Cold Block with Thermochromic PCR Cold Block Riser from -20°C freezer and thaw for **10 minutes** during Round 3 Plate Thaw.
- When prompted**, move the dual reservoir holder from the slanted plate holder back to Deck A.
- When prompted**, remove the 40 µm cell strainer.



**Note:** There may be bubbles left on the strainer. This will not affect the results.


- At the conclusion of the run, remove Parse Cold Block from the Deck C and place it on ice.

## 1.6. Round 3 Ligation


The program loads the Cell Suspension Mix on Deck A into the Round 3 plate on Deck B. Following the prompts for Round 3 incubation, the program will then load Round 3 Ligation Stop Buffer into all the wells and pool all the samples together on Deck A.

1. Gather the following components and reagents:


ITEM	SOURCE	QTY	HANDLING AND STORAGE
VIAFLO Pipette 12-Ch, 5-125 µL	INTEGRA Component	N/A	
Tip Deck for VIAFLO Pipette 12-Ch, 5-125 µL	INTEGRA Component	N/A	
Thermochromic PCR Cold Block	Parse-Provided	1	Pull the freezer block with riser from the -20°C freezer and leave them at room temperature for 10 minutes prior to use.
Thermochromic PCR Cold Block Riser	Parse-Provided	1	
INTEGRA Dual 25mL Basin Reservoir Adapter	INTEGRA Component	N/A	
10 mL transport tube	Consumables	1	
25 mL Basin Reservoir Liners	INTEGRA-Provided	2	
40 µm cell strainer	Consumables	1	
125 µL Tip Rack	INTEGRA-Provided	1	
1250 µL Tip Rack	INTEGRA-Provided	1	
○ Round 3 Stop Buffer	-20°C Reagents (Parse Reagents)	1	Thaw at room temperature then store on ice. Mix by vortexing.


2. Place Thermochromic PCR Cold Block with Thermochromic PCR Cold Block Riser on Deck B.
3. Replace the left reservoir.
4. Remove the  Round 3 Plate from the thermocycler and centrifuge for **1 minute** at 100 x g at 4°C.



- Place the  Round 3 Plate within the Thermochromic PCR Cold Block with Thermochromic PCR Cold Block Riser located on Deck B with A1 oriented towards the bottom left corner, and remove the seal. Deck layout should correspond to the configuration below.

## Deck Configuration

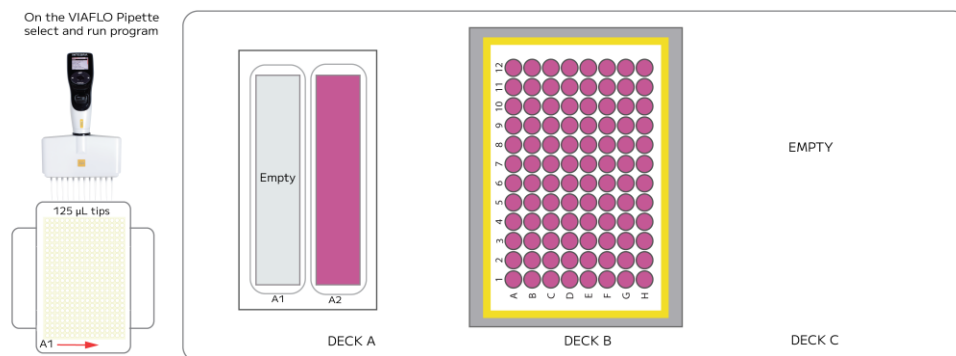

	DECK A	DECK B	DECK C
HARDWARE	Dual Reservoir Adapter	Thermochromic PCR Cold Block Thermochromic PCR Cold Block Riser	
CONSUMABLES	25 mL Reservoir Liner		
REAGENTS	● Samples	 Round 3 Plate	


- Remove D-ONE Pipetting Module 1-Ch, 5-1250  $\mu$ L and corresponding tip deck. Replace it with VIAFLO Pipette 12-Ch, 5-125  $\mu$ L and corresponding Tip Deck.



**Note:** Before removing the D-ONE Pipetting Module 1-Ch, 5-1250  $\mu$ L, ensure it is disconnected by pressing the back button until the main menu is displayed.

**Note:** Ensure the Tip Deck and Tip Boxes are securely positioned by pressing them down firmly until they click into place.

- Once the pipette and tip deck with tips are loaded on deck, select and run the program **WT S1 St6 V3** following the diagram below.

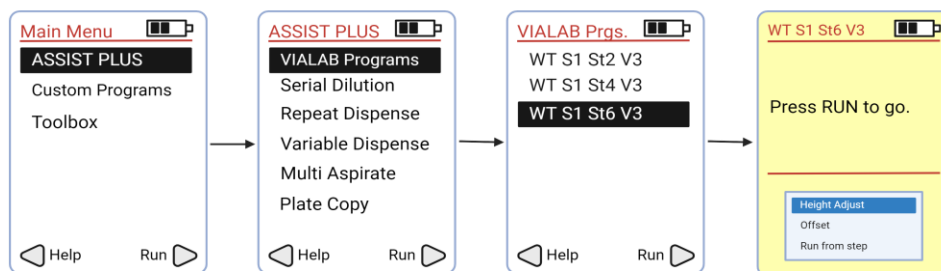

8. **When prompted**, reseal the ☐ Round 3 Plate with an adhesive seal and incubate the ☐ Round 3 Plate in a thermocycler with the following protocol.

BARCODING ROUND 3		
Run Time	Lid Temperature	Sample Volume
15 min	50°C	60 µL
Step	Time	Temperature
1	15 min	16°C
2	Hold	4°C

9. **When prompted**, replace the 25 mL basin reservoir liners with new 25 mL basin reservoir liners. Dispose used basin liners in biohazard waste.
10. **When prompted**, place the ☐ Round 3 Plate on Deck B,
11. **When prompted**, follow the prompts to add all the ☒ Round 3 Stop Buffer to the A2 basin ensuring that the volume added into the reservoir basin is evenly distributed.
12. Press "Run" to continue.
13. At the conclusion of the run, we recommend removing the Thermochromic PCR Cold Block with the Thermochromic PCR Cold Block Riser from the deck and storing them in the freezer.

## 1.7. Pre-Lysis


The pooled cell suspension on Deck B is strained into the 10 mL transport tube on Deck C4. The ● Spin Additive is then added into the cells and centrifuged. Supernatant is removed, the cells are resuspended in ○ Pre Lyse Wash Buffer, and centrifuged again. The supernatant is then removed.

1. Gather the following components and reagents:

ITEM	SOURCE	QTY	HANDLING AND STORAGE
D-ONE Pipetting Module 1-Ch, 5-1250 µL	INTEGRA Components	N/A	
Tip Deck for D-ONE Pipetting Module	INTEGRA Components	N/A	
Parse Cold Block	Parse-Provided	N/A	Keep on ice when not in use.
Dual 25mL Basin Reservoir Adapter	INTEGRA Components	N/A	
10 mL transport tube	Consumables	1	
25 mL basin reservoir liners	INTEGRA-Provided	2	
40 µm cell strainer	Consumables	1	
125 µL Tip Rack	INTEGRA-Provided	1	
1250 µL Tip Rack	INTEGRA-Provided	1	
● Spin Additive	4°C Reagents	1	Keep at room temperature.
○ Pre-Lysis Wash Buffer	-20°C Reagents	1	Thaw at room temperature then store on ice. Mix by pipetting 3x.

2. Add slanted plate holder on Deck B and incorporate a **10° tilt** (lowest side adjacent to Deck A).
3. Place the Parse Cold Block 1 on Deck C.
4. Place the reagent tubes in their respective orientation found in the deck configuration.
5. Put the 10 mL transport tube with the cell strainer in the C4 position of the Parse Cold Block 1. Deck layout should correspond to the configuration below.

## Deck Configuration

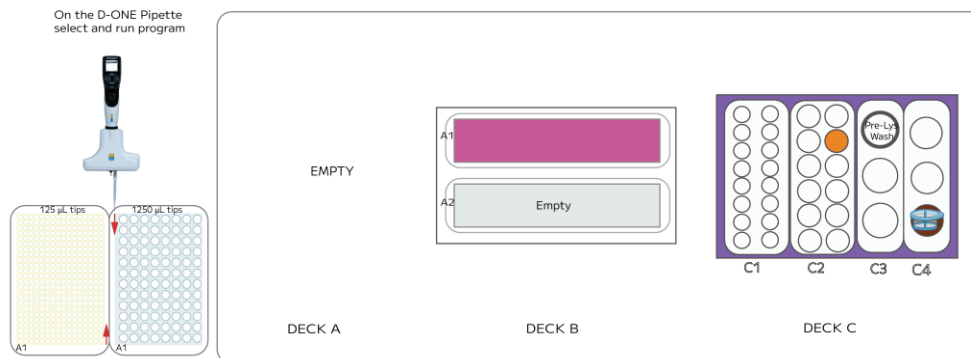

	DECK A	DECK B	DECK C
HARDWARE		Slanted Plate Holder (10°) Dual Reservoir Adapter	Parse Cold Block
CONSUMABLE		25 mL reservoir liner	10 mL transport tube and cell strainer
REAGENTS		● Samples	● Spin Additive ○ Pre Lyse Wash Buffer


- Remove the VIAFLO Pipette 12-Ch, 5-125 µL Pipette and corresponding Tip Deck. Attach the D-ONE Pipetting Module 1-Ch, 5-1250 µL and corresponding Tip Deck.



**Note:** Before removing the VIAFLO Pipette 12-Ch, 5-125 µL, ensure it is disconnected by pressing the back button until the main menu is displayed.

**Note:** Ensure the Tip Deck and Tip Boxes are securely positioned by pressing them down firmly until they click into place.

- Remove reagent caps, select and run the program **WT S1 St7 V3** following the diagram below.

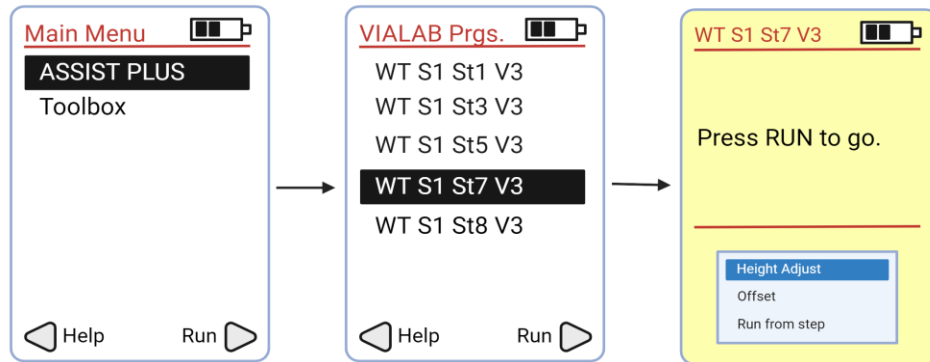

8. **When prompted**, move the dual reservoir holder back to Deck A.
9. **When prompted**, remove the cell strainer.
10. **When prompted**, cap and invert the 10 mL transport tube and centrifuge the pooled cells in a swinging bucket centrifuge cooled to 4°C for **10 minutes** at 250 x g.
11. **When prompted** after the centrifugation is complete, remove the cap from the 10 mL transport tube and place it back in its original position within the Parse Cold Block 1 located on Deck C4.
12. Press "Run" on Pipette.
13. **When prompted**, centrifuge the 10 mL transport tube in a swinging bucket centrifuge cooled to 4°C for **10 minutes** at 250 x g for a second spin.
14. **When prompted** after the centrifugation is complete, remove the cap from the 10 mL transport tube and place it back in its original position within the Parse Cold Block 1 located on Deck C4.
15. Press "Run" to continue.
16. While minimizing time on ice, mix and count the number of cells in the sample from the 10 mL transport tube on Deck C4 with a hemocytometer or alternative cell counting device. Record the cell count.
17. Clear the deck.

## 1.8. Lysis and Sublibrary Generation

The program requires a dilution of the sample to 500 cells per  $\mu\text{L}$  with a volume of 250  $\mu\text{L}$ . The program will create eight lysates with 12,500 cells each. It then mixes the Lysis Master Mix and aliquots the master mix into the generated lysates.

### To generate and lyse sublibraries:

1. Gather the following components and reagents:

ITEM	SOURCE	QTY	HANDLING AND STORAGE
D-ONE Pipetting Module 1-Ch, 5-1250 $\mu\text{L}$	INTEGRA Component	N/A	
Tip Deck for D-ONE Pipetting Module	INTEGRA Component	N/A	
Parse Cold Block	Parse-Provided	1	Keep on ice when not in use.
8 PCR strip tubes	Consumables	1	
1.5 mL tube	Consumables	2	
125 $\mu\text{L}$ Tip Rack	INTEGRA-Provided	1	
1250 $\mu\text{L}$ Tip Rack	INTEGRA-Provided	1	
● Pre-Lysis Dilution Buffer	-20°C Reagents	1	Thaw at room temperature then store on ice. Mix by pipetting 3x.
● Lysis Buffer	4°C Reagents	1	Place in a 37°C water bath until use.
● Lysis Enzyme	-20°C Reagents	1	Place directly on ice. Briefly centrifuge before use.

2. Dilute the cells to a concentration of 500 cells per  $\mu\text{L}$  for a total volume of **250  $\mu\text{L}$**  using the ● Pre-Lysis Dilution Buffer in a 1.5 mL transport tube. Store this dilution on ice until ready to use.




**Note:** The robot uses 25  $\mu\text{L}$  of this sample dilution to generate the lysates. If lysates of another size is required, dilute the sample suspension to a different concentration. (Example: A sample suspension with the concentration of 300 cells/nuclei per  $\mu\text{L}$  will create lysates with 7,500 cells/nuclei.)

3. Place a new PCR strip tube on Deck C1.


- Place the Lysis Buffer, Lysis Enzyme, empty 1.5 mL transport tube, and diluted sample from Section 1.8.2 on Deck C2. Deck should correspond to the Deck Configuration below.

### Deck Configuration

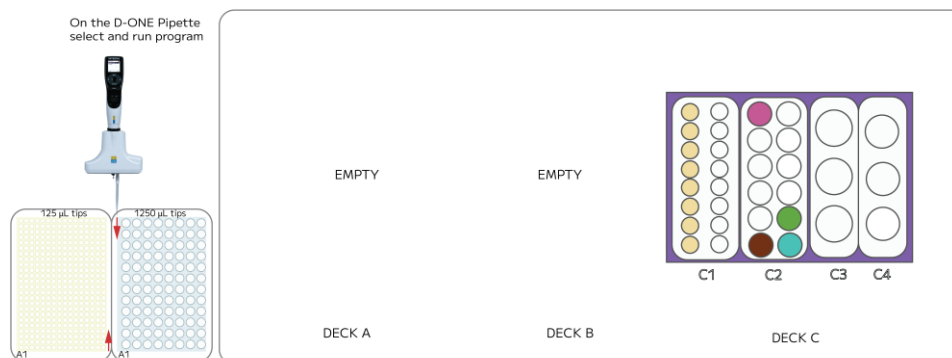

	DECK A	DECK B	DECK C
HARDWARE			Parse Cold Block
CONSUMABLES			<ul style="list-style-type: none"> <li>8-count PCR Tube</li> <li>1.5 mL Transport Tube</li> </ul>
REAGENTS			<ul style="list-style-type: none"> <li>Samples</li> <li>Lysis Buffer</li> <li>Lysis Enzyme</li> </ul>

- Remove the reagent caps, ensure that all strip tubes are open, select and run the program **WT S1 St8 V3** following the diagram below.

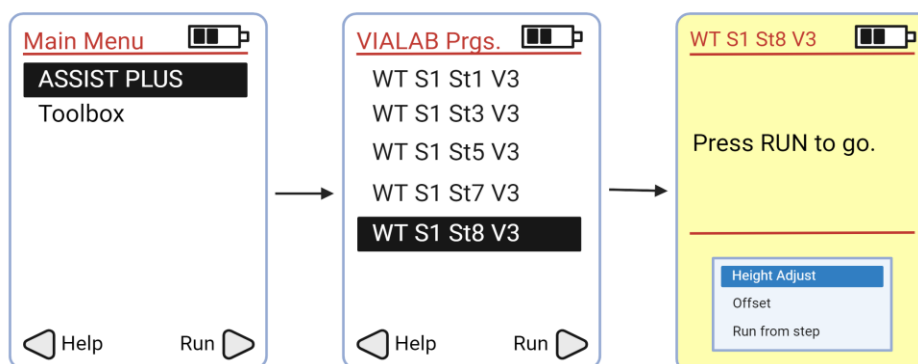

- Vortex the 0.2 mL tube(s) for **10 seconds**. Briefly centrifuge.

7. Place the tube(s) into a thermocycler and run the following program.

CELL LYSIS		
Run Time	Lid Temperature	Sample Volume
15 min	80°C	55 µL
Step	Time	Temperature
1	15 min	65°C
2	Hold	4°C

8. If continuing to Section 2 without freezing the sample, proceed to **Section 2: cDNA Capture and Amplification** in the [Evercode BCR User Guide](#) while the program is still running.



Safe stopping point: Sublibrary lysates can be stored at -80°C for up to 6 months.




## Appendices

### Appendix A: Pipetting Programs

#### Section 1.1. Sample Normalization


##### WT S1 St1 V3

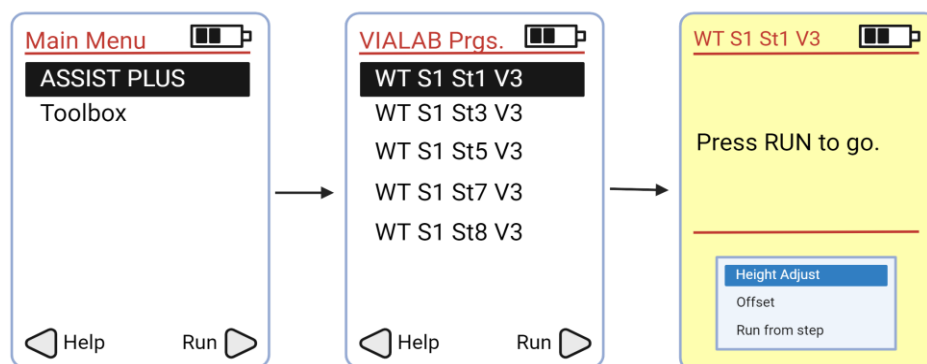

STEPS	ACTION
1	Initial Volumes
2	Diluent Worklist
3	Discard Tip
4	Sample Worklist
5	"Thaw R1 Plate" message

#### Section 1.2. Round 1 Plate Loading and Pooling


##### WT S1 St2 V3

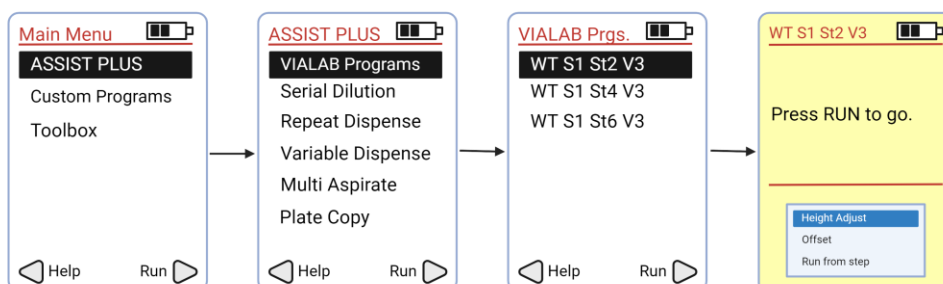

STEPS	ACTION
1	Initial Volumes
2	"Change Plate on Deck C to R1 plate" message
3	Transfer Diluted Fixed Sample to Round 1 Barcoding Plate
4	"Seal and incubate for Round 1 RT" message
5	"Plate R1 on Deck B" message
6	Volume Change
7	Pool Cells to Intermediate Wells

### Section 1.3. Round 2 Ligation Preparation

#### WT S1 St3 V3

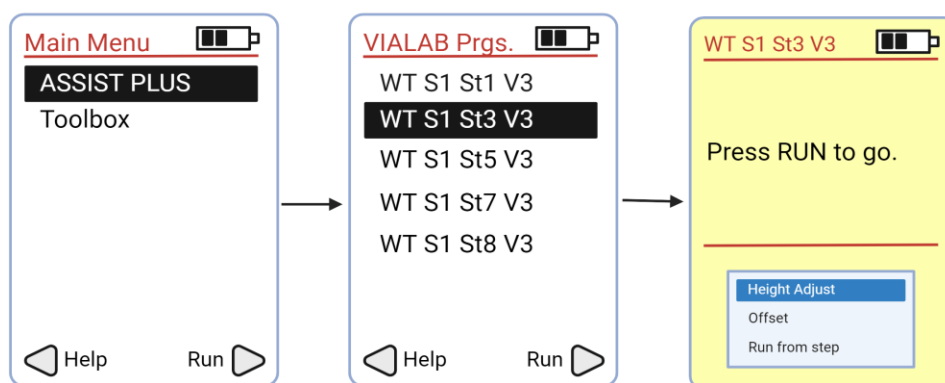


STEPS	ACTION
1	Initial Volumes
2	Pool row A into 10 mL Tube with 1000 $\mu$ L tips
3	Pool row into 10 mL Tube with 125 $\mu$ L tips
4	Add Spin Additive
5	"Invert Tube and Spin for 10 min" message
6	"Thaw R2 Plate" message

STEPS	ACTION
7	"Return tube to Deck C –" message
8-13	Remove Supernatant
14-15	Resuspend Cell Pellet with Resuspension Buffer
16	Add R2 Ligation Enzyme to Buffer
17	Mix Ligation Mix
18	Volume Change
19	Transfer Cells to Ligation Mix
20	Mix Cells in Ligation Mix
21-24	Transfer Ligation Mix to Basin
25-27	Mix Sample in Basin

## Section 1.4. Round 2 Ligation

### WT S1 St4 V3

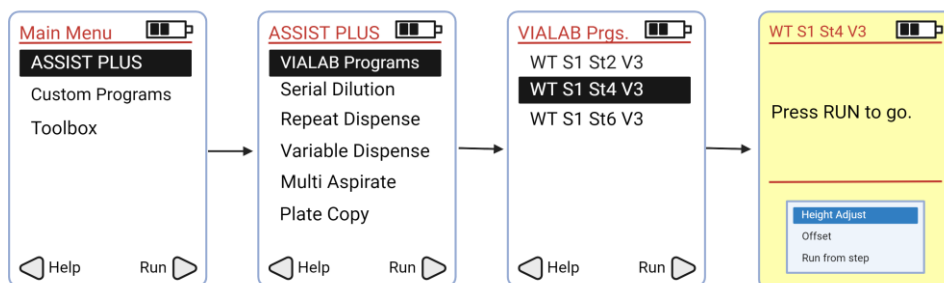


STEPS	ACTION
1	Initial Volumes
2-4	Mix samples in basins
5-12	Load Sample into Round 2 Plate
13	"Seal and incubate for R2 Ligation" message

STEPS	ACTION
14	"Replace both basin liners" message
15	"Reload R2 plate on Deck B" message
16	"Add R2 Stop to right basin" message
17	Volume Change
18	Add Stop
19	"Seal and incubate for Round 2 Stop" message
20	"Thaw R3 plate" message
21	"Replace right basin liner" message
22	"Reload R2 Plate on Deck B" message
23	Pool R2 Plate to basin

## Section 1.5. Round 3 Ligation Preparation

### WT S1 St5 V3

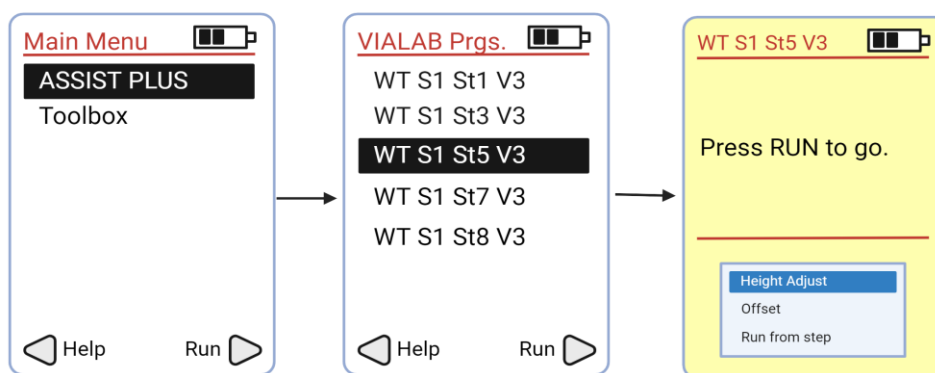


STEPS	ACTION
1	Initial Volumes

STEPS	ACTION
2	"Move basin to slanted holder" message
3	"Insert cell strainer in 10 mL tube" message
4	Wash basin mix
5	Volume change
6-10	Strain cells
11	"Move Basin Holder to Deck A" message
12	"Remove cell strainer" message
13	Volume change
14	Add R3 Ligation Enzyme
15	Mix Ligation Enzyme with sample
16	Volume change
17-21	Transfer cells to right basin

## Section 1.6. Round 3 Ligation

### WT S1 St6 V3

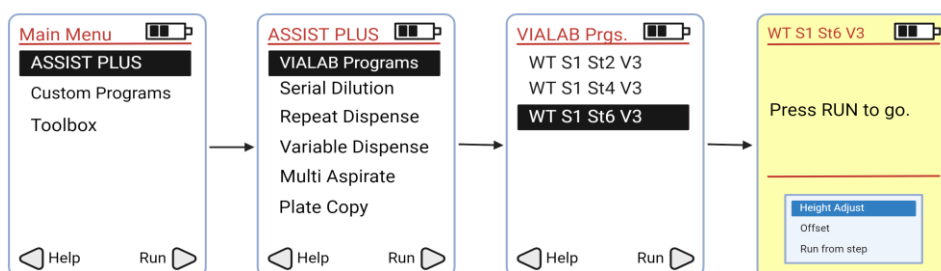


STEPS	ACTION
1	Initial Volumes

STEPS	ACTION
2-4	Mix sample in basin
5-12	Add sample to Round 3 Plate
13	"Seal and incubate for R3 Ligation" message
14	"Change both basin liners" message
15	"Reload R3 Plate on Deck B" message
16	"Add R3 Stop to right basin" message
17	Volume change
18	Add R3 Stop to plate
19	Pool R3 Plate

## Section 1.7. Pre-Lysis

### WT S1 St7 V3

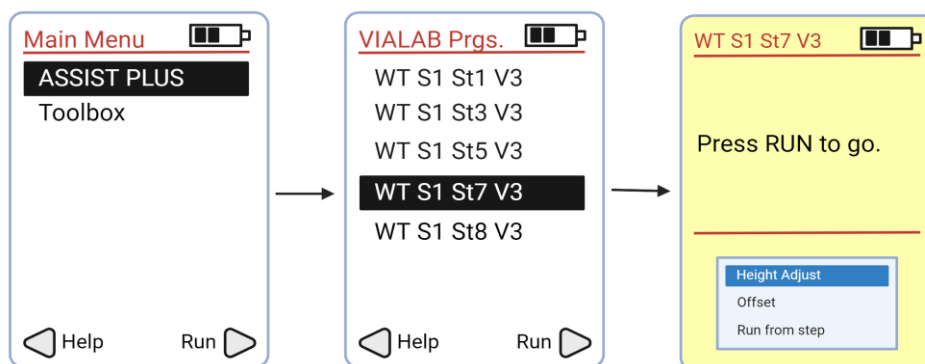


STEPS	ACTION
1	Initial Volumes
2	"Move Basin Holder to Deck B" message
3	"Insert cell strainer in 10 mL tube" message
4	Wash basin
5	Volume change
6-12	Strain cells

STEPS	ACTION
13	"Move Basin Holder to Deck A" message
14	"Remove cell strainer" message
15	Add Spin Additive
16	"Invert and spin for 10 min" message
17	"Return 10 mL tube to Deck C" message
18	Volume change
19-23	Remove supernatant
24-25	Resuspend Pre-Lyse
26	"Spin for 10 min" message
27	"Return 10 mL Tube to Deck C" message
28	Volume change
29-33	Remove supernatant
34	"Count cells/nuclei" message

## Section 1.8. Lysis and Sublibrary Generation

### WT S1 St8 V3

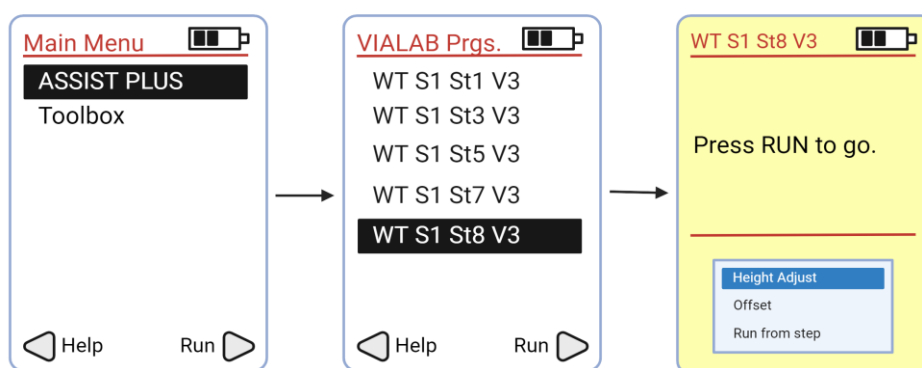

STEPS	ACTION
1	Initial Volumes

STEPS	ACTION
2-3	Aliquot 25 $\mu$ L of sample
4	Lysis Mastermix
5	Mix Lysis Mastermix
6	Add Lysis Mastermix to sample
7	"Vortex and centrifuge samples" message



## Appendix B: Troubleshooting

### Error warning during the execution of a program

- In case of an error warning, abort the program. Do not press "RUN". Scroll down "Run from step". You can then continue from the error location. This allows to resume directly from the point of error.

### Tip Pinching

- Tip pinching may occur when using a newly frozen PCR freezer block for a plate pooling step. Ensure that your PCR Freezer Block has been thawed at room temperature for at least 10 minutes.
- Ensure to apply sufficient pressure when placing the plate onto the frozen PCR freezer block.

### Incorrect tip location inputted

- The D-ONE Pipette will check the pipette tips before the step starts. The ASSIST PLUS will display an error warning if a pipette tip is incorrectly detected. Abort and restart the program with the correct tip location input.
- The VIAFLO will check the pipette tip after the step is finished. The ASSIST PLUS will display an error warning if a pipette tip is incorrectly detected. Abort and restart the program with the correct tip location input.

## Appendix C: Revision History

Version	Description	Date
1.0	Initial release	April 2024
1.1	Section 1.1: Updates Cell Number and Dilutions	May 2024
1.2	Updated configuration decks	June 2024
1.3	Updated Section 1	February 2025



[parsebiosciences.com](https://parsebiosciences.com)

[support@parsebiosciences.com](mailto:support@parsebiosciences.com)

